



INTERVENTIONAL  
VASCULAR  
DIAGNOSTICS  
AND THERAPY

## Surecan<sup>®</sup> Safety II

Power-injectable safety needles for Access Ports

# Surecan® Safety II needles for Access Ports



## Surecan® Safety II is a non-coring safety needle with a small size and unique design.

The Surecan® Safety II guarantees easy handling and maximal safety to minimize needlestick injuries. This ensures a high comfort for clinicians and patients, in hospitals as well as for home care treatment.

### Indicated for up to seven days for:

- Infusion of chemotherapy
- Total parenteral nutrition (TPN)
- Antibiotic therapy
- Transfusion
- Pain relief and drainage
- Blood sampling
- Power injections of contrast media up to 325 psi

### Ergonomic low profile wings

- Provide secure handling during puncture and removal
- Permit small and comfortable dressing
- Needle gauge indicated on wings

### Closed-cell pad

- Increases patient comfort

### Special needle tip

- Aids needle retention in the septum during power injections

### Unique bevel

- Non-coring bevel designed to minimize septum damage

### Coloured clamp

- Possible interruption of infusion at any time
- Easy to identify gauge of the needle
- Contains maximum psi information for high pressure injections

### Extension tubing

- Transparent anti-kink tubing
- Power injectable

Extensive choice of references:

19 G to 22 G

12 mm to 38 mm



Surecan® Safety II, compatible with all B. Braun Access Port Systems

## Surecan® Safety II allows hospitals to reply to the requirements of the European Council Directive 2010/32/EU from 10th May 2010: "Prevention from sharp injuries"

- "Health and safety of workers is paramount and is closely linked to the health of patients."\*
- Surecan Safety II is a "medical (device) ... incorporating safety-engineered protection mechanisms."\*


\* European council directive 2010/32/EU from 10th May 2010

### Green point and "click" confirms safety mechanism

- Easy confirmation of correct deployment
- Double affirmation through visual and audible feature




### Transparent base

- Easier visualisation of the insertion point
- Holes aid skin aeration

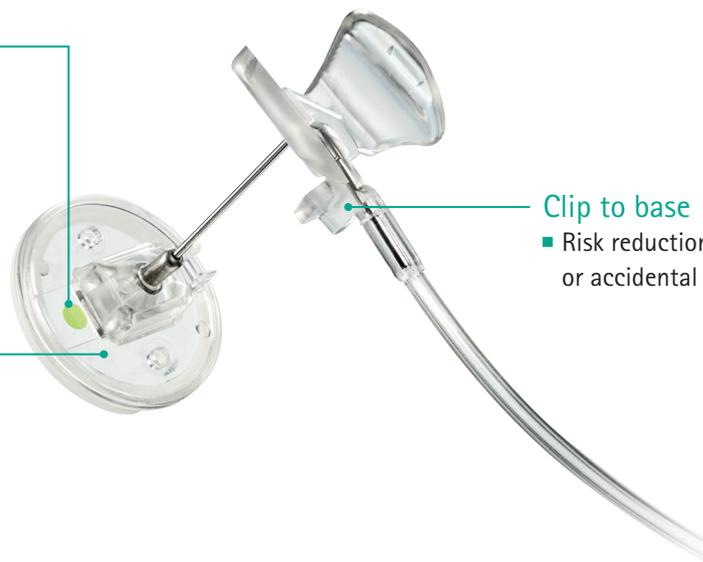

### Clip to base

- Risk reduction of separation or accidental detachment

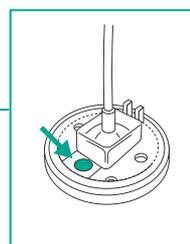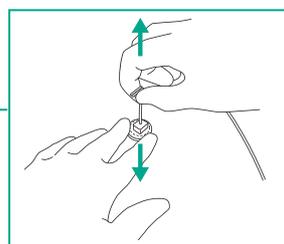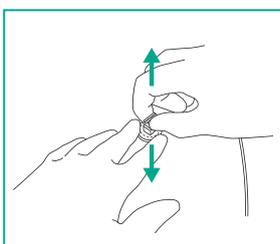

### Safe removal


- Stabilise the base of the needle
- Firmly pull the wings up until you hear a „click“ and see the green point

For more information on the risks and prevention of needle-stick injuries, please visit [www.safeinfusiontherapy.com](http://www.safeinfusiontherapy.com) and refer to our brochure no. 6069095.

## Caresite® as pre-connected needle-free valve


The Surecan® Safety II is also available with the pre-connected needle-free valve Caresite®. The split septum valve technology combined with the ergonomic design enables clear and easy access for improved patient care.

- Protects against microbiological contamination
- Prevents catheter occlusion through positive displacement (valve design)
- Ease of use

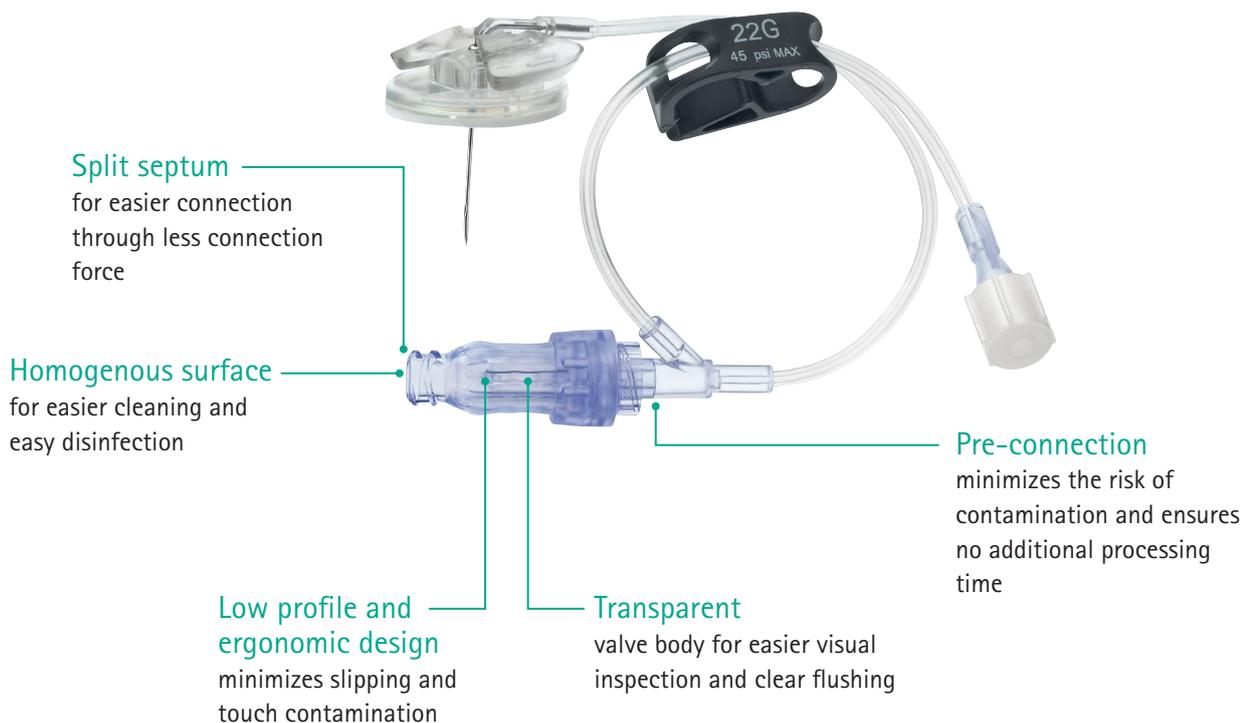

### ATTENTION

High pressure (HP) injections have to be performed through the Caresite® valve

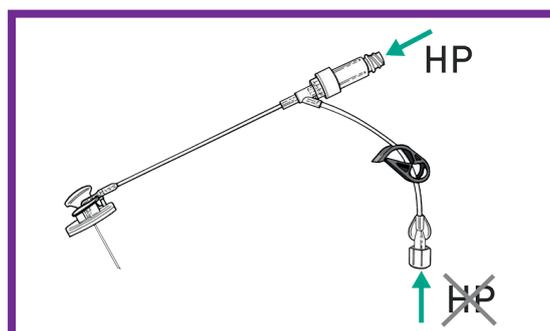

## Surecan® Safety II safety needle

### Surecan® Safety II

|   |  G19 | REF                        |  G20 | REF      |  G22 | REF      | Unit       |
|---|---|----------------------------|---|----------|---|----------|------------|
| Priming volume/ml   | 0.30  |                            | 0.25  |          | 0.20  |          |            |
|  | 12 mm   | 04447042                   | 12 mm   | 04447043 | 12 mm   | 04447044 | 20 per box |
|   | 15 mm   | 04447000                   | 15 mm   | 04447005 | 15 mm   | 04447010 | 20 per box |
|   | 20 mm   | 04447001                   | 20 mm   | 04447006 | 20 mm   | 04447011 | 20 per box |
|   | 25 mm   | 04447002                   | 25 mm   | 04447007 | 25 mm   | 04447012 | 20 per box |
|   | 32 mm   | 04447003                   | 32 mm   | 04447008 | 32 mm   | 04447013 | 20 per box |
|   | 38 mm   | 04447004                   | 38 mm   | 04447009 | -   | -        | 20 per box |
| <b>DEHP</b><br>FREE   | <b>LATEX</b><br>FREE  | <b>POWER</b><br>INJECTABLE | <b>MR</b><br>CONDITIONAL  |          |   |          |            |

### Surecan® Safety II with Y-site and Caresite®

|  |  G19 | REF                        |  G20 | REF      |  G22 | REF      | Unit       |
|--|---|----------------------------|---|----------|---|----------|------------|
| Priming volume/ml  | 0.65  |                            | 0.55  |          | 0.45  |          |            |
|  | 12 mm   | 04447057                   | 12 mm   | 04447058 | 12 mm   | 04447059 | 20 per box |
|  | 15 mm   | 04447045                   | 15 mm   | 04447050 | 15 mm   | 04447054 | 20 per box |
|  | 20 mm   | 04447046                   | 20 mm   | 04447051 | 20 mm   | 04447055 | 20 per box |
|  | 25 mm   | 04447047                   | 25 mm   | 04447052 | 25 mm   | 04447056 | 20 per box |
|  | 32 mm   | 04447048                   | 32 mm   | 04447053 | -   | -        | 20 per box |
|  | 38 mm   | 04447049                   | -   | -        | -   | -        | 20 per box |
| <b>DEHP</b><br>FREE  | <b>LATEX</b><br>FREE  | <b>POWER</b><br>INJECTABLE | <b>MR</b><br>CONDITIONAL  |          |   |          |            |

## Accessories

### Caresite®

needle free connector for high pressure injections

|   | REF       | Description     | Unit        |
|---|-----------|-----------------|-------------|
|  | 415122-01 | Caresite® valve | 100 per box |

### SwabCap®

small safety device for disinfection of needle-free membrane valves

|   | REF      | Description                           | Unit        |
|---|----------|---------------------------------------|-------------|
|  | EM-SCXT3 | SwabCap® disinfection device 70 %-IPA | 200 per box |

### Omniflush® with SwabCap®

0,9 % NaCl prefilled syringe for flushing with disinfection device (only for application in the non-sterile area)

|   | REF          | Description                             | Unit        |
|---|--------------|---|-------------|
|  | EM-3513572SC | Omniflush® with SwabCap® 3 ml in 10 ml  | 100 per box |
|   | EM-3513575SC | Omniflush® with SwabCap® 5 ml in 10 ml  | 100 per box |
|   | EM-35135726C | Omniflush® with SwabCap® 10 ml in 10 ml | 100 per box |

### Medibox®

medical sharps collector for easy and safe disposal

|   | REF     | Description    | Unit |
|---|---------|----------------|------|
|  | 9193405 | Medibox® 0.7 l | 1    |
|   | 9193413 | Medibox® 3.0 l | 1    |
|   | 9193421 | Medibox® 5.0 l | 1    |

### Needlebox

for secure storage of non-safety items

|   | REF     | Description                      | Unit |
|---|---------|----------------------------------|------|
|  | 8910670 | Needlebox for non-safety needles | 1    |

Distributor

B. Braun Melsungen AG | Vascular Systems | Sieversufer 8 | 12359 Berlin | Germany  
Phone +49 30 568207-300 | Fax +49 30 568207-210 | [www.bbraun.com](http://www.bbraun.com)

Manufacturer acc. to MDD 93/42/EEC

B. Braun Medical | 204 Avenue du Maréchal Juin | 92100 Boulogne-Billancourt | France

The product trademarks 'Surecan', 'Ultrasite', 'Caresite', 'SwabCap', 'Omniflush' and 'Medibox' are registered trademarks of B. Braun Melsungen AG.

Subject to technical changes. All rights reserved. This brochure may only be used for the exclusive purpose of obtaining information about our products. Reproduction in any form partial or otherwise is not permitted.