



AESCULAP® Sterile Container System

General Catalog

AESCULAP® Sterile Container System

Content

Sterile Containers

- PRIMELine Pro
- Basis
- Mini sterile containers
- Dental sterile containers
- Sterile containers for scopes
- XXL endo sterile container and accessories
- Da Vinci SP Camera Rack
- Disposal containers

Instrument Baskets

- Basket configurator
- Basket with feet
- Basket without feet
- Basket inserts
- Lid for baskets JF, JG
- Special trays
- Scope trays
- Storage rack for MIC instruments and trocars

Aesculap InOS – Instrument Organization System

- Instrument racks and holders
- Silicon mats
- Storage solution for uterine dilators
- Mounting instructions
- Holding clamps
- Fixation elements
- Stop frames

Accessories

- Identification labels
- Filter accessories

Consumables

- Reusable filters
- Paper filters
- Indicator labels
- Seals
- Maintenance products
- Labels, bar code labels
- Labels, printer accessories

Often imitated. Never equaled

The AESCULAP® original



Always better. From the beginning.

In 1971, Aesculap revolutionized the packaging and supply of sterile instruments with the world's first sterile container. From then up to the present day, the Aesculap original has become the universal solution and Aesculap is the global sterile container market leader.



AESCULAP® Sterile Container System

Leading global supplier for decades

Aesculap, a leading manufacturer of sterile containers for over 50 years.



With the development and launch of the first real sterile container in 1971, Aesculap delivered a world first. Their innovative technology and their ease of handling have made AESCULAP® sterile containers a permanent fixture in hospitals worldwide. The use of AESCULAP® sterile containers means that the entire process can be standardized, controlled and documented: from storage to transport, through to the aseptic presentation of the sterile instruments in the operating room. Even at the disposal stage, the AESCULAP® sterile container system reliably protects users in health facilities from contamination.


The AESCULAP® sterile container system provides customers with an extensive and clinically proven product portfolio for sterile packaging.

Various sterile container models and sizes, as well as an extensive range of accessories, comprising different basket models, together with the proven Aesculap InOS® organization system enable reliable protection for almost all surgical instrument sets.





The right instruments
at the right time
in sterile, functional condition
at the right place of use
at a high level of profitability¹
in a reliable packaging system

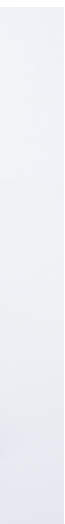

¹ Krohn M, Fengler J, Mickley T, Flessa S. Analysis of processes and costs of alternative packaging options of sterile goods in hospitals - a case study in two German hospitals. Health Economics Review (2019) 9:1

AESCULAP® Sterile Container System

Leading global supplier for decades



Perforated soft packaging can have serious consequences ...



... focus more on stability with the AESCULAP® Sterile Container System.

¹ Krohn M, Fengler J, Mickley T, Flessa S. Analysis of processes and costs of alternative packaging options of sterile goods in hospitals – a case study in two German hospitals. Health Economics Review (2019) 9:1

² Waked WR, Simpson AK, Miller CP, Magit DP, Grauer JN. Sterilization wrap inspections do not adequately evaluate instrument sterility. Clin Orthop Relat Res. 2007; 462:207–211

AESCULAP® Sterile Container System

Quality without compromise – Made in Germany

Innovation Factory in Tuttlingen – Innovation meets tradition

The AESCULAP® sterile container system has been manufactured at the Tuttlingen/Germany site for more than 50 years. All the components of an AESCULAP® sterile container comply with the national and international requirements and standards (e.g. DIN EN ISO 11607-1). The solid, robust design* and the manufacturing technology used by Aesculap enable the sterile containers, instrument baskets and all the other components to be manufactured with very precise dimensions and a high level of stability.





With a modern industrial architecture, at the Innovation Factory, which was opened in June 2015, Aesculap acts on the B. Braun corporate values: be innovative, increase efficiency and live sustainably. The goal is an intelligent factory. This is characterized by adaptability, efficient use of resources and ergonomic design.

Think about tomorrow today!

The Innovation Factory meets the criteria for the German Sustainable Building Council (DGNB) silver certificate. Among other things, the factory has a fully automated

container production line with 10 processing stations, which are linked with seven robots, as well as a fully automated raw materials warehouse with approx. 600 storage areas.

The spacious premises with the latest equipment has space for around 500 employees. Production, development and product management are spatially linked to one another. This ensures a quick and constant exchange between all departments – the basis for a future-oriented finding and implementation of solutions.

* EN 868-8 (handle strength test, stacking test)

AESCULAP® Sterile Container System

Quality without compromise – Made in Germany

AESCULAP® sterile containers are tested based on the latest national and international standards.



EN ISO 11607-1
EN 868-8



Aesculap Technical Service



The AESCULAP® sterile containers protect sterile items against external influences and keep them sterile for when they are used. If they are damaged in any way (e.g. through deformation or wear), they can usually be repaired.

Original AESCULAP® container spare parts are compatible with all generations so the sterile containers will always be upgraded to the latest status quo by a repair. As a result, every repair by the Aesculap Technical Service helps retain the value of your sterile container stock.

Don't use damaged sterile containers! Please check your AESCULAP® sterile container before each use. Precise instructions can be found in the Aesculap instructions for use TA015110 on <https://eifu.bbraun.com> and on the Aesculap functional testing posters C63302, C92802 (please ask us about them if needed).

AESCULAP® Sterile Container System

High flexibility: As the bottoms are identical, you can choose between two different systems. Using existing AESCULAP® container bottoms helps reduce purchasing costs.

The extensive and well-proven Aesculap product portfolio for sterile instrument packaging offers a range of packaging solutions for all hospital instrument sets used in practice. Sterile containers in various models and sizes and the comprehensive range of accessories make the Sterile Container System a universal solution for the sterile supply of surgical instruments.



PRIMELine® PRO



BASIS




AESCULAP® Sterile Container System

PRIMELine® PRO

The PRIMELine® PRO aluminum container lid is equipped with Aesculap's proven permanent germ barrier system. This system consists of a cartridge with a well-protected, integrated PTFE (polytetrafluoroethylene) filter. This germ barrier has been validated for 5,000 sterilization cycles according to EN 868-8 (2009), Annex H and therefore commonly provides lifelong function (in everyday hospital use).

As a result, typically, no filter changes are necessary, thereby saving time and process costs. A cleaning validation has been carried out for the PRIMELine® Pro lid in an automated procedure. A positive cleaning result, suitable for everyday use in the hospital, was achieved*.

The cleaning result fully complied with the acceptance criteria.

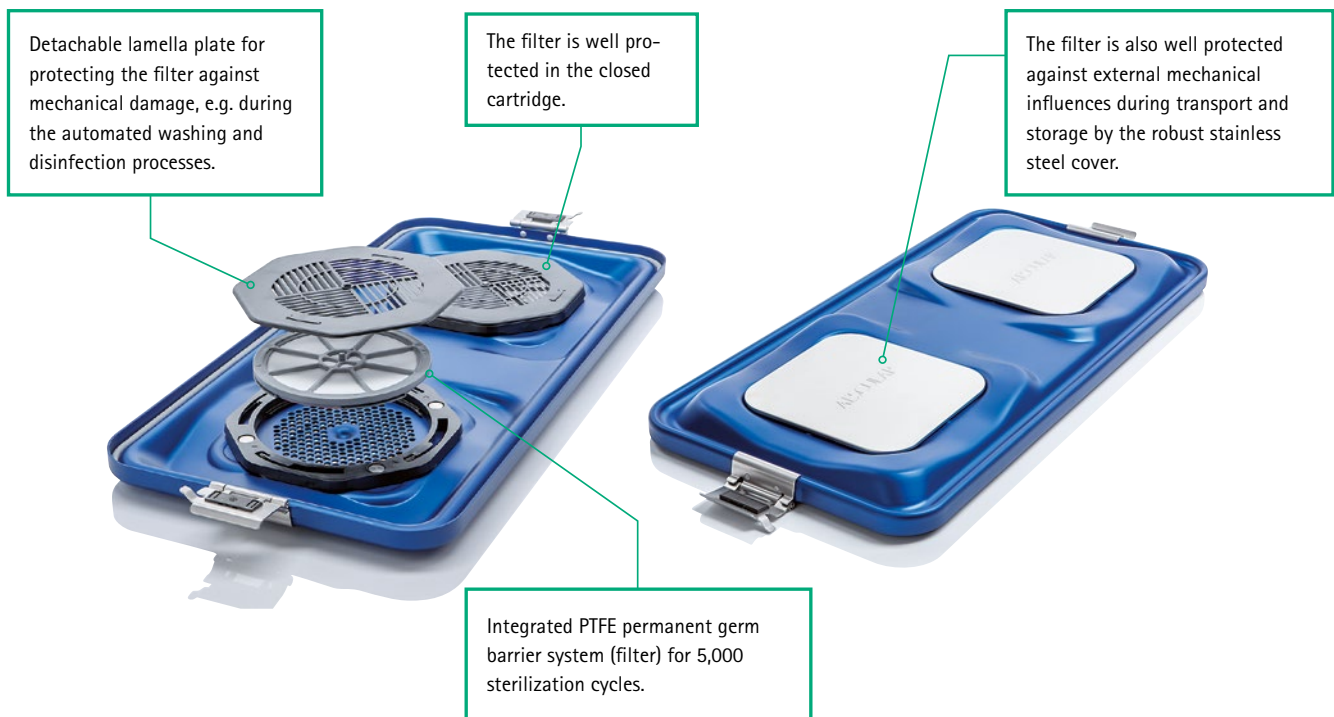

- Stable, robust aluminum lid
- Typically no filter change – saves time and process costs
- Cleaning validation fully complies with the acceptance criteria
- Available in 5 colors (red, blue, green, gold and silver)
- Compatible with all container bottoms (1/1, 1/2, 3/4) manufactured since 1988

* Please refer to the brochure PRIMELine® Pro Test Results C01702

Mode of action of an AESCULAP® sterile container in steam sterilization. For more information please refer to D-ST22016. Please ask us about it if needed.



BASIS

The basic model made of aluminum for use with AESCULAP® single-use sterile filters or with AESCULAP® PTFE reusable filters for 1,000 sterilizations.

The Basis version can be fitted with a protective lid, if desired. However, the sterile container is still fully functional without a protective lid. The two lids can only be separated when the sterile container is opened for better cleaning. There is a push-button on the inside of the lid for this purpose.



AESCLAP® Sterile Container System

AESCLAP® Sterile Filters

Different filters for different requirements



Aesculap PTFE (polytetrafluorethylene) permanent germ barrier system validated for 5,000 sterilization cycles according to EN 868-8, Annex H for lifelong function (in everyday hospital use). As a result, typically, no filter changes are necessary, saving time and process costs. For steam sterilization in a validated fractionated vacuum process according to EN ISO17665-1. Integrated into the PRIMELine® PRO sterile container lid.

AESCLAP® paper filter for single use. The filters for steam sterilization have an indicator according to EN ISO 11140-1 Type 1. Therefore, used and unused filters can be clearly differentiated visually. For use in Basis sterile container lids.

AESCLAP® reusable PTFE (polytetrafluorethylene) filters for 1,000 sterilization cycles (for steam sterilization in a fractionated vacuum process validated according to EN ISO17665-1). For use in Basis sterile container lids.



Mode of action of filters used for sterile containers within steam sterilization.
For more information please refer to D-ST20144. Please ask us about it if needed.



Paper filters and reusable filters are held firmly in the lid by the AESCULAP® filter retention plate (Basis). It has a silicone gasket and ensures that the filter is securely fitted. The retention plate has a safety lock (two release buttons) to prevent an uncontrolled release of the filter. It also enables quick filter changes. The audible click makes control easier.



Filter retention plate with safety lock


- Prevents uncontrolled release
- Quick filter change
- Audible click-in



AESCULAP® Sterile Container System


The AESCULAP® sterile container system in daily use In the medical device reprocessing unit:

After cleaning and disinfection, a visual function check is performed and, if necessary, a new paper filter inserted.

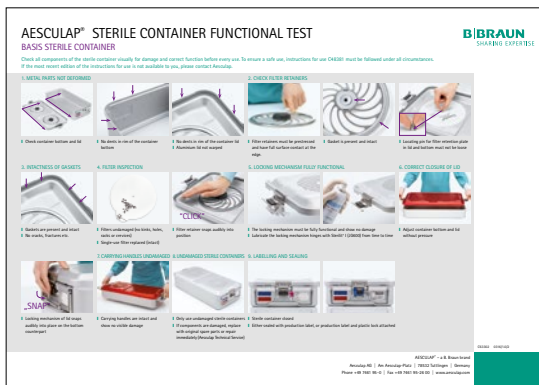

Next, the sterile container is then sealed on both sides. The production label is then affixed. The sterile container is now ready for sterilization.

The sterile container can be labeled using the color-coded labels, e.g. with information on the contents of the basket and/or where it belongs (neuro, OR 1, building A, etc.).

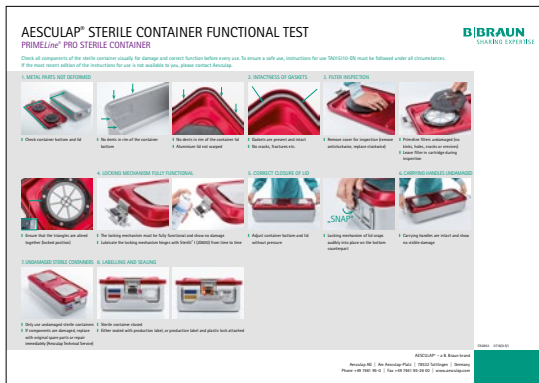

The packed instrument basket is then placed in the container bottom. The lid is then put on and locked.

Sterile container sealed with the AESCULAP® plastic seal with indicator¹.



Sterile container sealed with AESCULAP® indicator label¹.

¹ According to EN ISO 11140-1 Type 1

In the operating room

After sterilization and the corresponding cooling phase, the sterile containers are placed in the sterile instrument store. The instrument sets are now ready to be used again. The AESCULAP® sterile container storage period (up to 1 year) has been tested in several long-term studies. The tests proved that sterility was maintained throughout the entire period. The storage conditions used in the tests complied with ANSI/AAMI ST79².

The corresponding nursing staff bring the closed sterile containers to the operating room and check that the seals are intact and the indicator has changed. The sterile container is then opened. When it is opened, the container seals are destroyed by the locking flaps*. A nurse lifts the lid off and stands back. A check is then carried out to ensure that the filter and gasket are intact. The lid is then put down.

The surgical nurse picks the instrument set up in the sterile container and places it on the preparation table.



* Alternatively:
Break the seal
by hand.

² ANSI/AAMI ST79: Comprehensive Guide to Steam Sterilization and Sterility Assurance in Health Care Facilities

AESCULAP® Sterile Container System

Quick, reliable and efficient

Which sterile container is suitable for each sterilisation process? – An overview
For all AESCULAP® sterile containers with CE marking, listed in this catalog.

		Steam sterilization ¹		Low temperature processes		
		Fractionated vacuum process	Gravity process	Sterrad® 100S, 200, NX, 100NX	Steris® Amsco V-Pro 1, V-Pro 1 Plus	Eto ²
Sterile containers	Basis	JK Series ✓	JN Series ✓	Incompatible	Incompatible	Not validated
	Vario	JK Series ✓	JN Series ✓	Incompatible	Incompatible	Not validated
	PRIMELine® PRO	JK Series ✓	Not validated	Incompatible	Incompatible	Incompatible
	Scope containers	✓	Not validated	JM020+JM021 ✓	JM020+JM021 ✓	✓
	Mini containers	JK Series ✓	JN Series ✓	JM188+JM174 ✓	JM188+JM174 ✓	Not validated
	Dental containers	JK Series ✓	JN Series ³ ✓	Incompatible	Incompatible	Not validated
	XXL EndoContainer	JK Series ✓	JN Series ✓	Incompatible	Incompatible	Not validated
Filters	Paper – for single use	✓	✓	Incompatible	Incompatible	✓
	PTFE – reusable	✓	✓	Incompatible	Incompatible	Not validated
	Polypropylene – for single use (JF167)	Not validated	Not validated	✓	✓	Not validated

OK (sometimes only with certain article numbers – see note)
Not validated
Incompatible

¹ The sterilization must be done in a steam sterilization process validated according to ISO 17665-1 in a sterilizer according to EN 285.

² The sterilization must be validated according to EN 550/ISO 11135-1.

³ Except JN292, JN294, JN295, JN296 – these are for the fractionated vacuum process only.

Aesculap Academy

Connect. Exchange. Enable.



We drive lifelong learning for health care professionals.

The Aesculap Academy, as the continuing education forum of the B. Braun Group, offers scientific and product-specific training and workshops in the healthcare sector at its locations in Tuttlingen and Berlin, Germany, as well as in 43 other countries worldwide – both virtually and in person, in accordance with our motto: “Connect. Exchange. Enable.”




Vision – Mission – Commitment

Our Vision: We set sustainable standards for the benefit of patients and to preserve our environment for future generations.

Our Mission: The AESCULAP® Sterile Container System provides the right instrument, in a sterile condition and with low process costs, at the right time in the operating room.

Our Commitment: We have been in dialog with our customers for more than 50 years in order to provide them with reliable solutions. This dialog will be at the heart of our business in the future as well.

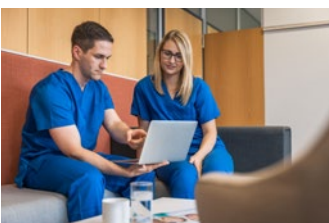

AESCLAP[®] Sterile Container System

Stepping up to the sustainability challenge in healthcare

Healthcare professionals around the world are looking for new ways to make their activities and the medical service they provide more sustainable and environmentally friendly. As the climate challenge comes into ever sharper focus, hospitals are endeavoring to do their part to improve sustainability and environmental performance in the healthcare sector.

For further information about
AESCLAP Aicon[®]

Scan the QR code or visit our Website:
www.bbraun.com/aesclapaicon





Sterile Containers

PRIMELine® Pro

Size 1/1

Lids for individual color combination

	Color
JP101	red
JP102	blue
JP103	green
JP104	gold
JP105	silver



	
	JK440
	JK441
	JK442
	JK444
	JK446

Solid bottoms = JK

Tables of dimensions PRIMELine® Pro

Size 1/1

Bottom 1/1		Lid		Sterile container with lid
	Outside dimensions		Outside dimensions	 Outside dimensions Inside dimensions
JK440	592 x 274 x 90	JP101 - JP105	588 x 281 x 36	592 x 281 x 107 544 x 258 x 75
JK441	592 x 274 x 120			592 x 281 x 137 544 x 258 x 105
JK442	592 x 274 x 135			592 x 281 x 152 544 x 258 x 120
JK444	592 x 274 x 187			592 x 281 x 204 544 x 258 x 172
JK446	592 x 274 x 247			592 x 281 x 264 544 x 258 x 232

PRIMELine® Pro

Size 3/4

Lids for individual color combination

	Color
JP111	red
JP112	blue
JP113	green
JP114	gold
JP115	silver



	
	JK740
	JK741
	JK742
	JK744

Solid Bottoms = JK

Tables of dimensions PRIMELine® Pro

Size 3/4

Bottom 3/4		Lid		Sterile container with lid
	Outside dim.		Outside dim.	 Outside dimensions Inside dimensions
JK740	470 x 274 x 90	JP111 - JP115	465 x 281 x 36	470 x 281 x 107 421 x 258 x 75
JK741	470 x 274 x 120			470 x 281 x 137 421 x 258 x 105
JK742	470 x 274 x 135			470 x 281 x 152 421 x 258 x 120
JK744	470 x 274 x 187			470 x 281 x 204 421 x 258 x 172

PRIMELine® Pro

Size 1/2

Lids for individual color combination

 for JK339	 for JK340-JK346	Color
JP131	JP121	red
JP132	JP122	blue
JP133	JP123	green
JP134	JP124	gold
JP135	JP125	silver



	
	JK339 without carrying handle
	JK340
	JK341
	JK342
	JK344
	JK346

Solid bottoms = JK

Tables of dimensions PRIMELine® Pro

Size 1/2

Bottom 1/2		Lid		Sterile container with lid
	Outside dimensions		Outside dimensions	Outside dimensions Inside dimensions
JK339 without carrying handle	300 x 274 x 57	JP131 - JP135	298 x 281 x 36	300 x 281 x 74 260 x 265 x 35
JK340	300 x 274 x 90	JP121 - JP125	298 x 281 x 36	300 x 281 x 107 254 x 258 x 75
JK341	300 x 274 x 120			300 x 281 x 137 254 x 258 x 105
JK342	300 x 274 x 135			300 x 281 x 152 254 x 258 x 120
JK344	300 x 274 x 187			300 x 281 x 204 254 x 258 x 172
JK346	300 x 274 x 247			300 x 281 x 264 254 x 258 x 232

Basis

Size 1/1

Lids for individual color combination

Inner lid

	Color
JK485	red
JK486	blue
JK487	green
JK488	gold
JK489	silver

Outer lid

	Color
JK475A	red
JK476A	blue
JK477A	green
JK478A	gold
JK479A	silver



		
	JK440	JN440
	JK441	JN441
	JK442	JN442
	JK444	JN444
	JK446	JN446

Solid bottoms = JK

Perforated bottoms = JN



JF101R support for trays

Tables of dimensions Basis

Size 1/1

Bottom 1/1		Inner lid		Outer lid		Sterile container with inner lid		Sterile container with inner and outer lid	
	Outside dim.		Outside dim.		Outside dim.		Outside dim. Inside dim.		Outside dim. Inside dim.
JK440	592 x 274 x 90	JK485 - JK489	588 x 285 x 36	JK475A - JK479A	588 x 291 x 36	592 x 285 x 108	592 x 291 x 115	544 x 258 x 75	544 x 258 x 75
JK441	592 x 274 x 120					592 x 285 x 138	592 x 291 x 145	544 x 258 x 105	544 x 258 x 105
JK442	592 x 274 x 135					592 x 285 x 153	592 x 291 x 160	544 x 258 x 120	544 x 258 x 120
JK444	592 x 274 x 187					592 x 285 x 205	592 x 291 x 212	544 x 258 x 172	544 x 258 x 172
JK446	592 x 274 x 247					592 x 285 x 265	592 x 291 x 272	544 x 258 x 232	544 x 258 x 232

Bottom 1/1		Inner lid		Outer lid		Sterile container with inner lid		Sterile container with inner and outer lid	
	Outside dim.		Outside dim.		Outside dim.		Outside dim. Inside dim.		Outside dim. Inside dim.
JN440	592 x 274 x 90	JK485 - JK489	588 x 285 x 36	JK475A - JK479A	588 x 291 x 36	592 x 285 x 108	592 x 291 x 115	544 x 258 x 63	544 x 258 x 63
JN441	592 x 274 x 120					592 x 285 x 138	592 x 291 x 145	544 x 258 x 93	544 x 258 x 93
JN442	592 x 274 x 135					592 x 285 x 153	592 x 291 x 160	544 x 258 x 108	544 x 258 x 108
JN444	592 x 274 x 187					592 x 285 x 205	592 x 291 x 212	544 x 258 x 160	544 x 258 x 160
JN446	592 x 274 x 247					592 x 285 x 265	592 x 291 x 272	544 x 258 x 220	544 x 258 x 220

Basis

Size 3/4


Lids for individual color combination

Inner lid

	Color
JK785	red
JK786	blue
JK787	green
JK788	gold
JK789	silver

Outer lid

	Color
JK775A	red
JK776A	blue
JK777A	green
JK778A	gold
JK779A	silver

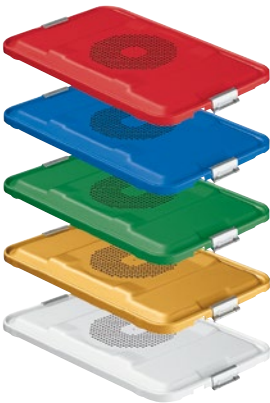

		
	JK740	JN740
	JK741	JN741
	JK742	JN742
	JK744	JN744

Solid bottoms = JK

Perforated bottoms = JN



JF104R support for trays

Tables of dimensions Basis

Size 3/4

Bottom 3/4		Inner lid		Outer lid		Sterile container with inner lid		Sterile container with inner and outer lid	
	Outside dim.		Outside dim.		Outside dim.		Outside dim. Inside dim.		Outside dim. Inside dim.
JK740	470 x 274 x 90	JK785 - JK789	460 x 285 x 36	JK775A - JK779A	460 x 291 x 36	470 x 285 x 108	421 x 258 x 75	470 x 291 x 115	421 x 258 x 75
JK741	470 x 274 x 120					470 x 285 x 138	421 x 258 x 105	470 x 291 x 145	421 x 258 x 105
JK742	470 x 274 x 135					470 x 285 x 153	421 x 258 x 120	470 x 291 x 160	421 x 258 x 120
JK744	470 x 274 x 187					470 x 285 x 205	421 x 258 x 172	470 x 291 x 212	421 x 258 x 172

Bottom 3/4		Inner lid		Outer lid		Sterile container with inner lid		Sterile container with inner and outer lid	
	Outside dim.		Outside dim.		Outside dim.		Outside dim. Inside dim.		Outside dim. Inside dim.
JN740	470 x 274 x 90	JK785 - JK789	460 x 285 x 36	JK775A - JK779A	460 x 291 x 36	470 x 285 x 108	421 x 258 x 63	470 x 291 x 115	421 x 258 x 63
JN741	470 x 274 x 120					470 x 285 x 138	421 x 258 x 93	470 x 291 x 145	421 x 258 x 93
JN742	470 x 274 x 135					470 x 285 x 153	421 x 258 x 108	470 x 291 x 160	421 x 258 x 108
JN744	470 x 274 x 187					470 x 285 x 205	421 x 258 x 160	470 x 291 x 212	421 x 258 x 160

Basis

Size 1/2

Lids for individual color combination

Inner lid

Outer lid

 for JK339	 for JK340-JK346	Color	 for JK339-JK346	Color			
JK395	JK385	red	JK375A	red		JK339 without carrying handle	-
JK396	JK386	blue	JK376A	blue		JK340	JN340
JK397	JK387	green	JK377A	green		JK341	JN341
JK398	JK388	gold	JK378A	gold		JK342	JN342
JK399	JK389	silver	JK379A	silver		JK344	JN344
						JK346	JN346



Solid bottoms = JK

Perforated bottoms = JN



JF102R support for trays

Tables of dimensions Basis

Size 1/2

Bottom 1/2		Inner lid		Outer lid		Sterile container with inner lid	Sterile container with inner and outer lid
	Outside dim.		Outside dim.		Outside dim.	 Outside dim. Inside dim.	 Outside dim. Inside dim.
JK339	300 x 274 x 57 without carrying handle	JK395 - JK399	298 x 285 x 36	JK375A - JK379A	298 x 291 x 36	300 x 285 x 67 260 x 265 x 35	300 x 291 x 75 260 x 265 x 35
JK340	300 x 274 x 90	JK385 - JK389	298 x 285 x 36	JK375A - JK379A	298 x 291 x 36	300 x 285 x 108 254 x 258 x 75	300 x 291 x 115 254 x 258 x 75
JK341	300 x 274 x 120					300 x 285 x 138 254 x 258 x 105	300 x 291 x 145 254 x 258 x 105
JK342	300 x 274 x 135					300 x 285 x 153 254 x 258 x 120	300 x 291 x 160 254 x 258 x 120
JK344	300 x 274 x 187					300 x 285 x 205 254 x 258 x 172	300 x 291 x 212 254 x 258 x 172
JK346	300 x 274 x 247					300 x 285 x 265 254 x 258 x 232	300 x 291 x 272 254 x 258 x 232

Bottom 1/2		Inner lid		Outer lid		Sterile container with inner lid	Sterile container with inner and outer lid
	Outside dim.		Outside dim.		Outside dim.	 Outside dim. Inside dim.	 Outside dim. Inside dim.
JN340	300 x 274 x 90	JK385 - JK389	298 x 285 x 36	JK375A - JK379A	298 x 291 x 36	300 x 285 x 108 254 x 258 x 63	300 x 291 x 115 254 x 258 x 63
JN341	300 x 274 x 120					300 x 285 x 138 254 x 258 x 93	300 x 291 x 145 254 x 258 x 93
JN342	300 x 274 x 135					300 x 285 x 153 254 x 258 x 108	300 x 291 x 160 254 x 258 x 108
JN344	300 x 274 x 187					300 x 285 x 205 254 x 258 x 160	300 x 291 x 212 254 x 258 x 160
JN346	300 x 274 x 247					300 x 285 x 265 254 x 258 x 220	300 x 291 x 272 254 x 258 x 220

Mini Sterile Containers

Mini sterile containers for small sets and small instruments

The mini sterile containers were developed for small sets (e.g. small-wound treatment) and small, delicate instruments especially for microsurgery.

Suitable sterilisation processes (please also refer to page 20):

- Mini sterile containers are suitable for steam sterilisation using the fractionated vacuum process (solid bottom) and using the gravity process (perforated bottom).
- Separately labeled mini sterile containers (please also refer to the next page) are additionally suitable for low-temperature sterilisation processes with hydrogen peroxide in a Sterrad® 100 S, Sterrad® 200, Sterrad® NX, Sterrad® 100 NX, Steris V-Pro® 1 and Steris V-Pro® 1 Plus.
- Mini sterile containers are not suitable for hot air sterilization!



Mini Sterile Containers

Solid bottom = JK



Perforated bottom = JN



Lids for individual color combination



		Outside dimensions Length x Width x Height (mm)
JK187	JN187	310 x 132 x 30
JK188	JN188	310 x 132 x 57
	JM188*	310 x 132 x 57

	Color	Outside dimensions Length x Width x Height (mm)
JK170	red	305 x 140 x 25
JK171	blue	305 x 140 x 25
JK172	green	305 x 140 x 25
JK173	gold	305 x 140 x 25
JK174	silver	305 x 140 x 25
JM174*	silver	305 x 140 x 25

Dimensions complete - Length x Width x Height (mm)

	Sterile containers, complete  Outside dimensions Inside dimensions
Bottom	with Lid
JK187	310 x 140 x 40 274 x 126 x 21
JK188	310 x 140 x 67 274 x 126 x 48

	Sterile containers, complete  Outside dimensions Inside dimensions
Bottom	with Lid
JN187	310 x 140 x 40 274 x 126 x 18
JN188	310 x 140 x 67 274 x 126 x 45
JM188*	310 x 140 x 67 274 x 126 x 45

* For low-temperature sterilisation processes with hydrogen peroxide in a Sterrad® 100 S, Sterrad® 200, Sterrad® NX, Sterrad® 100 NX, Steris V-Pro® 1 and Steris V-Pro® 1 Plus. To be used with disposable filter JF167 (see page 101).

Dental Sterile Containers

Dental Sterile Containers – the design not only for dental surgery

AESCULAP® dental sterile containers are particularly suitable for sterile preparation in dental surgeries. The Dental Sterile Container fits into common small sterilizers. A wide range of accessories is available for combining individual instrument sets.



- The well-established locking levers facilitate the aseptic presentation
- Sterile containers can be stacked for space-saving in the sterile supplies cycle
- Quick-release filter retention plate for quick and easy replacement of AESCULAP® sterile filters
- System-designed trays for processing, thermal disinfection, sterilization, transport and storage of instruments
- Suitable sterilisation processes (please also refer to page 20): Dental sterile containers are suitable for steam sterilisation using the fractionated vacuum process (solid bottom) and using the gravity process (perforated bottom). Dental sterile containers are not suitable for hot air sterilization!

Note | This sterile container can be ordered complete (lid and body) under one article number.

Solid bottom

	Outside dim. Inside dim.		Outside dim.		Color	Outside dim.
	JN292 312 x 190 x 46 277 x 176 x 27	JK286	312 x 185 x 37	JN091	silver	319 x 190 x 25
	JN294 312 x 190 x 65 277 x 176 x 46	JK288	312 x 185 x 56			
	JN295 312 x 190 x 92 277 x 176 x 73	JK289	312 x 185 x 83			
	JN296 312 x 190 x 130 277 x 176 x 111	JK290	312 x 185 x 121			

Perforated bottom

	Outside dim. Inside dim.		Outside dim.		Color	Outside dim.
	JN092 312 x 190 x 46 277 x 176 x 20	JN086	312 x 185 x 37	JN091	silver	319 x 190 x 25
	JN094 312 x 190 x 65 277 x 176 x 39	JN088	312 x 185 x 56			
	JN095 312 x 190 x 92 277 x 176 x 66	JN089	312 x 185 x 83			
	JN096 312 x 190 x 130 277 x 176 x 104	JN090	312 x 185 x 121			

Sterile Containers for Scopes

Separate storage and provision of the scopes

- Separate, quick aseptic supply, independent of instrument sets
- Protection of sensitive scopes including light conducting cable during reprocessing, transportation and storage
- Individual design with the organisation system Aesculap InOS®
- Protected perforation field
- The handle on the tray lid allows simple aseptic removal
- Stackable
- Suitable sterilisation processes (please also refer to page 20):
The scope sterile container is suitable for steam sterilisation using the fractionated vacuum process. Separately labeled scope sterile containers (please also refer to the next page) are additionally suitable for low-temperature sterilisation processes with hydrogen peroxide in a Sterrad® 100 S, Sterrad® 200, Sterrad® NX, Sterrad® 100 NX, Steris V-Pro® 1 and Steris V-Pro® 1 Plus. The scope container is not suitable for using in gravity processes!



Note

Sterrad® 100 S, Sterrad® 200, Sterrad® NX and Sterrad® 100 NX are registered trademarks of ASP.
V-Pro® 1 and V-Pro® 1 Plus are registered trademarks of Steris.

Scope Sterile Containers

Lid

For the fractionated vacuum process
(incl. filter retention plate JK098)



	Outside dimensions Length x Width x Height (mm)
JK020	540 x 152 x 28

AESCULAP® reusable filter for steam
sterilisation, for use with filter retention
plate JK098, service life 1,000 cleaning
and sterilisation cycles.

Order unit (PAK): Pack of 2



	Outside dimensions Length x Width (mm)
JK091	95 x 215

Tray

For up to two laparoscopic scopes and one
light cable.
Preconfigured trays see page 61



	Outside dimensions Length x Width x Height (mm)
JF435R	501 x 134 x 45

Bottom

For the fractionated vacuum process
(incl. filter retention plate JK098)



	Outside dimensions Length x Width x Height (mm)
JN021	540 x 144 x 65


AESCULAP® disposable filter for steam
sterilisation with process indicator,
for use with filter retention plate JK098.

Order unit (PAK): Pack of 100



	Outside dimensions Length x Width (mm)
JF166	118 x 235

Identification label with adhesive back and
process indicator* for steam sterilisation and
sterilisation with ethylene oxide.
Order unit (PAK): Pack of 100

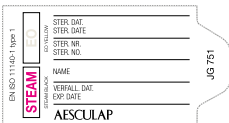

	Outside dimensions Length x Width (mm)
JG751	67 x 35

Sterile container, complete

Outside dimensions
Inside dimensions

540 x 152 x 84
503 x 134 x 45

Lid

For the hydrogen peroxide sterilisation in Sterrad® 100 S, Sterrad® 200, Sterrad® NX, Sterrad® 100 NX, Steris V-Pro® 1 and Steris V-Pro® 1 Plus, (incl. filter retention plate JK098)



Outside dimensions Length x Width x Height (mm)
JM020 540 x 152 x 28

Bottom

For the hydrogen peroxide sterilisation in Sterrad® 100 S, Sterrad® 200, Sterrad® NX, Sterrad® 100 NX, Steris V-Pro® 1 and Steris V-Pro® 1 Plus, (incl. filter retention plate JK098)



Outside dimensions Length x Width x Height (mm)
JM021 540 x 144 x 65

Sterile container, complete

Outside dimensions
Inside dimensions

540 x 152 x **84**
503 x 134 x **45**

Disposable filter for the hydrogen peroxide sterilisation in Sterrad® 100 S, Sterrad® 200, Sterrad® NX, Sterrad® 100 NX, Steris V-Pro® 1 and Steris V-Pro® 1 Plus.

Order unit (PAK): Pack of 100



Outside dimensions Length x Width (mm)
JF167 235 x 118


Universal filter retention plate



Outside dimensions Length x Width x Height (mm)
JK098 245 x 100 x 6

Identification label with adhesive back and process indicator* for vaporized hydrogen peroxide sterilisation.

Order unit (PAK): Pack of 100

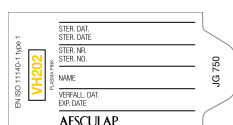

Outside dimensions Length x Width (mm)
JG750 67 x 35

XXL Endo Sterile Container and Accessories

- Extra long Sterile Container for endoscopic instruments with sheath length of up to 420 mm.
- The instruments can be sterilized as fully assembled units, if applicable.
- Suitable sterilisation processes (please also refer to page 20):
The XXL Endo Sterile Container is suitable for steam sterilisation using the fractionated vacuum process (solid bottom) and using the gravity process (perforated bottom).
- **Caution:** The Endo Sterile Container can be sterilized only in sterilizers of an appropriate chamber size (at least 720 mm effective depth).



Bottom

solid bottom = JK

perforated bottom = JN

		Outside dimensions Length x Width x Height (mm)		Outside dimensions Length x Width x Height (mm)
JK443	JN443	711 x 274 x 120	JK490	708 x 280 x 35
	JN445	711 x 274 x 187		

Lid

Dimensions complete -

Length x Width x Height (mm)


	Sterile container, complete
	Outside dimensions Inside dimensions
JK443	711 x 285 x 138 667 x 258 x 105
JN443	711 x 285 x 138 667 x 258 x 93
JN445	711 x 285 x 205 667 x 258 x 160

Baskets



	Outside dim. Inside dim.
JF232R	660 x 253 x 56 658 x 251 x 42
JF236R	660 x 253 x 86 658 x 251 x 72


Silicone positioning mat for storing of accessories

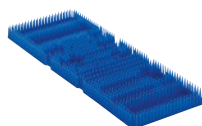

	Outside dim.
JF933	610 x 230 x 30

Endorack fitting for JN445

Storage rack for MIS-instruments, XXL comprising: storage rack, XXL, perforated tray and silicone mat for instruments with shaft length up to 420 mm, fitting for JN445 + JK490

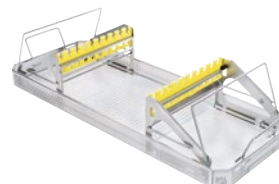


	Outside dim.
PL961R	660 x 253 x 166

Da Vinci SP* Camera Rack



Sterile Containers

The Da Vinci SP* Camera rack is used with the Endocontainer JN445 bottom and JK490 lid.

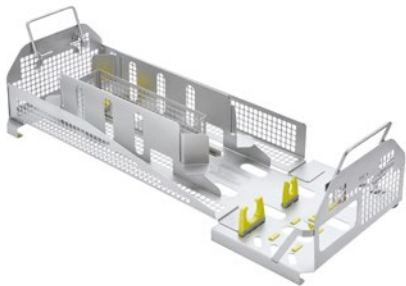

	Outside dim.
	Inside dim.
JF449R	660 x 253 x 163

*Registered trademark by Intuitive Surgical

Disposal Containers

Disposal containers complete the instrument cycle.
The closed containers protect both staff and the environment from contamination during the transport of used instruments.

The disposal containers are of the same dimensions as the standard sterile containers, thus allowing standardized and automated instrument return and transport.

As all components of the disposal containers are made of stainless steel, combined cleaning and disinfecting agents can be used in them, and they can be cleaned with alkaline detergents.



Container Size	Solid bottom		Lid with opening		Lid without opening		Container, complete  Outside dim. Inside dim.
		Outside dim.		Outside dim.		Outside dim.	
1/2	JK060R	295 x 275 x 90	JK067R	295 x 280 x 36	JK068R	295 x 280 x 36	295 x 280 x 102 260 x 260 x 80
1/2	JK062R	295 x 275 x 135	JK067R	295 x 280 x 36	JK068R	295 x 280 x 36	295 x 280 x 150 260 x 260 x 125
1/1	JK262R	586 x 275 x 135	JK267R	586 x 280 x 36	JK268R	586 x 280 x 36	586 x 280 x 150 550 x 260 x 125



Instrument Baskets

Basket Configurator

Illustrative overview for JF baskets

Which basket goes with which container bottom?



JF series

- Black PEEK plastic feet for stacking
- Fold-down handles on the side
- Mat surface
- 4 mm base perforation / 2 mm side perforation
- Can be combined with Aesculap InOS®
- Sloped edges

1/1 short


			Top basket								
			Basket in combination with:								
			Alone	JF211R	JF211R with lid*	JF212R	JF212R with lid*	JF213R	JF213R with lid*	JF214R	JF214R with lid*
Bottom basket	JF211R	Container bottom	JK440 JN440	JK441 JN441	JK441 JN441	JK441 JN442	JK442 JN442	JK442 JN444	JK444 JN444	JK444 JN444	JK444 JN444
	JF211R with lid*				JK441 JN442	JK442 JN444	JK442 JN444			JK444 JN444	JK444 JN446
	JF212R		JK440 JN441	JK441 JN442	JK442 JN442	JK442 JN444	JK444 JN444	JK444 JN444	JK444 JN446	JK444 JN446	JK446 JN446
	JF212R with lid*		JK441 JN441	JK442 JN444	JK442 JN444						
	JF213R		JK441 JN442		JK444 JN444	JK444 JN444	JK444 JN444	JK444 JN446	JK446 JN446	JK446 JN446	
	JF213R with lid*		JK442 JN444	JK444 JN444							JK444 JN444
	JF214R		JK444 JN444	JK444 JN444		JK444 JN446	JK446 JN446	JK446 JN446	JK446 JN446	JK446	
	JF214R with lid*			JK444 JN446		JK446 JN446	JK446 JN446		JK446	-	

*JF217R lid



1 Example

- Basket alone: JF212R (without lid)
- ➡ suitable container bottom: JK440 or JN441

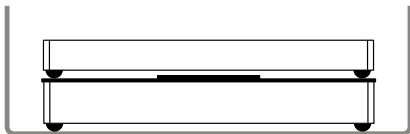

2 Example

- Top basket: JF212R (without lid)
- Bottom basket: JF213R with lid
- ➡ suitable container bottom: JK444 or JN444

JK = Solid bottom
(fractionated vacuum process)



JN = Perforated bottom
(gravity process)



Note | According to EN ISO 868-8, the load must not exceed a total of 10 kg in order to avoid excessive condensate formation.
Caution! The inner height and therefore the appropriate basket changes depending on whether JK or JN bottoms are used.



JF series

- Black PEEK plastic feet for stacking
- Fold-down handles on the side
- Mat surface
- 4 mm base perforation / 2 mm side perforation
- Can be combined with Aesculap InOS®
- Sloped edges


1/1 long

Top basket

Basket in combination with:


		Alone	JF221R	JF221R with lid*	JF222R	JF222R with lid*	JF223R	JF223R with lid*	JF224R	JF224R with lid*
Bottom basket	JF221R	JK440 JN440	JK441 JN441	JK441 JN441	JK441 JN442	JK442 JN442	JK442 JN444	JK444 JN444	JK444 JN444	JK444 JN444
	JF221R with lid*									
	JF222R	JK440 JN441	JK441 JN442	JK442 JN442	JK442 JN444	JK444 JN444	JK444 JN444	JK444 JN444	JK444 JN446	JK444 JN446
	JF222R with lid*									
	JF223R	JK441 JN442	JK442 JN444	JK442 JN444	JK444 JN444	JK444 JN444	JK444 JN446	JK444 JN446	JK446 JN446	JK446 JN446
	JF223R with lid*									
	JF224R	JK444 JN444	JK444 JN446	JK444 JN444	JK444 JN446	JK446 JN446	JK446 JN446	JK446 JN446	JK446	JK446
	JF224R with lid*									

*JF227R lid

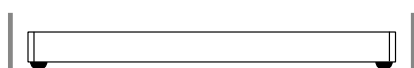

1 Example

- Basket alone: JF222R (without lid)
- ➔ suitable container bottom: JK440 or JN441

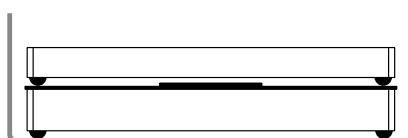

2 Example

- Top basket: JF222R (without lid)
- Bottom basket: JF223R with lid
- ➔ suitable container bottom: JK444 or JN444

JK = Solid bottom
(fractionated vacuum process)



JN = Perforated bottom
(gravity process)



Note | According to EN ISO 868-8, the load must not exceed a total of 10 kg in order to avoid excessive condensate formation.

Caution! The inner height and therefore the appropriate basket changes depending on whether JK or JN bottoms are used.

Basket Configurator

Illustrative overview for JF baskets

Which basket goes with which container bottom?




JF series

- Black PEEK plastic feet for stacking
- Fold-down handles on the side
- Mat surface
- 4 mm base perforation / 2 mm side perforation
- Can be combined with Aesculap InOS®
- Sloped edges

3/4


Top basket											
Basket in combination with:											
			Alone	JF251R	JF251R with lid*	JF252R	JF252R with lid*	JF253R	JF253R with lid*	JF254R	JF254R with lid*
Bottom basket	JF251R	Container bottom	JK740 JN740	JK741 JN741	JK741 JN741	JK741 JN742	JK742 JN742	JK742 JN744		JK744 JN744	JK744 JN744
	JF251R with lid*				JK741 JN742	JK742 JN744	JK742 JN744		JK744 JN744	JK744	JK744
	JF252R		JK740 JN741	JK741 JN742	JK742 JN742	JK742 JN744		JK744 JN744			-
	JF252R with lid*		JK741 JN741		JK742 JN744		JK744 JN744		JK744	-	-
	JF253R		JK741 JN742	JK742 JN744		JK744 JN744			JK744	-	-
	JF253R with lid*		JK742 JN744	JK744 JN744	JK744 JN744			JK744	-	-	-
	JF254R					JK744	-	-	-	-	-
	JF254R with lid*		JK744 JN744	JK744		-	-	-	-	-	-

*JF257R lid

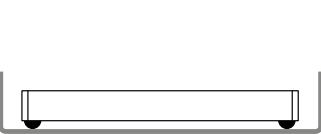

1 Example

- Basket alone: JF252R (without lid)
- ➔ suitable container bottom: JK740 or JN741

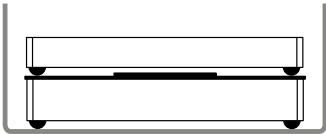

2 Example

- Top basket: JF252R (without lid)
- Bottom basket: JF253R with lid
- ➔ suitable container bottom: JK744 or JN744

JK = Solid bottom
(fractionated vacuum process)



JN = Perforated bottom
(gravity process)



Note | According to EN ISO 868-8, the load must not exceed a total of 7.5 kg in order to avoid excessive condensate formation.
Caution! The inner height and therefore the appropriate basket changes depending on whether JK or JN bottoms are used.



JF series

- Black PEEK plastic feet for stacking
- Fold-down handles on the side
- Mat surface
- 4 mm base perforation / 2 mm side perforation
- Can be combined with Aesculap InOS®
- Sloped edges


1/2

Top basket

Basket in combination with:


		Alone	JF111R	JF111R with lid*	JF112R	JF112R with lid*	JF113R	JF113R with lid*	JF114R	JF114R with lid*
Bottom basket	JF111R	Container bottom	JK340 JN340	JK341 JN341	JK341 JN341	JK341 JN342	JK342 JN342	JK342 JN344	JK344 JN344	JK344 JN344
	JF111R with lid*		JK340 JN340	JK341 JN341	JK341 JN342	JK342 JN342	JK342 JN344	JK342 JN344	JK344 JN344	JK344 JN346
	JF112R		JK340 JN341	JK341 JN342	JK342 JN342	JK342 JN344	JK344 JN344	JK344 JN344	JK344 JN346	JK346 JN346
	JF112R with lid*		JK341 JN341	JK342 JN344	JK342 JN344	JK344 JN344			JK346 JN346	
	JF113R		JK341 JN342	JK344 JN344	JK344 JN344	JK344 JN344	JK344 JN346	JK346 JN346	JK346 JN346	JK346 JN346
	JF113R with lid*		JK342 JN344	JK344 JN344						
	JF114R		JK344 JN344	JK344 JN346	JK344 JN346	JK346 JN346	JK346 JN346	JK346 JN346	JK346 JN346	JK346 JN346
	JF114R with lid*		JK344 JN346	JK344 JN346	JK346 JN346	JK346 JN346	JK346 JN346	JK346 JN346	JK346 JN346	JK346 JN346

*JF117R lid

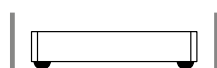

1 Example

- Basket alone: JF112R (without lid)
- ➔ suitable container bottom:
JK340 or JN341

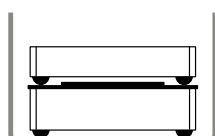

2 Example

- Top basket: JF112R (without lid)
- Bottom basket: JF113R with lid
- ➔ suitable container bottom:
JK344 or JN344

JK = **Solid** bottom
(fractionated vacuum process)



JN = **Perforated** bottom
(gravity process)




Note | According to EN ISO 868-8, the load must not exceed a total of 5 kg in order to avoid excessive condensate formation.

Caution! The inner height and therefore the appropriate basket changes depending on whether JK or JN bottoms are used.

Baskets

Basket **with** Feet

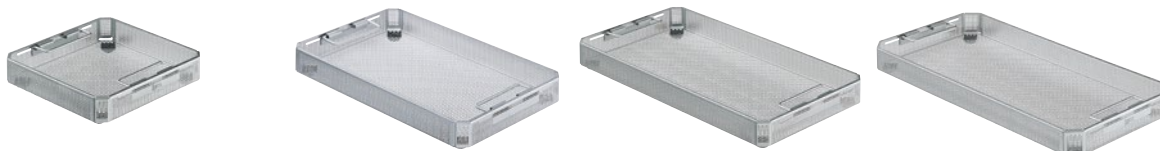

	 Standard 1/2	Outside dim. Inside dim.	 Standard 3/4	Outside dim. Inside dim.	 Standard 1/1 short	Outside dim. Inside dim.	 Standard 1/1 long	Outside dim. Inside dim.
 36	JF111R	243 x 253 x 36 241 x 251 x 22	JF251R	406 x 253 x 36 404 x 251 x 22	JF211R	485 x 253 x 36 483 x 251 x 22	JF221R	540 x 253 x 36 538 x 251 x 22
 56	JF112R	243 x 253 x 56 241 x 251 x 42	JF252R	406 x 253 x 56 404 x 251 x 42	JF212R	485 x 253 x 56 483 x 251 x 42	JF222R	540 x 253 x 56 538 x 251 x 42
 76	JF113R	243 x 253 x 76 241 x 251 x 62	JF253R	406 x 253 x 76 404 x 251 x 62	JF213R	485 x 253 x 76 483 x 251 x 62	JF223R	540 x 253 x 76 538 x 251 x 62
 106	JF114R	243 x 253 x 106 241 x 251 x 92	JF254R	406 x 253 x 106 404 x 251 x 92	JF214R	485 x 253 x 106 483 x 251 x 92	JF224R	540 x 253 x 106 538 x 251 x 92

Basket **without** Feet



	 Standard 1/2	Outside dim. Inside dim.	 Standard 3/4	Outside dim. Inside dim.	 Standard 1/1 short	Outside dim. Inside dim.	 Standard 1/1 long	Outside dim. Inside dim.
 24	JG111R	243 x 253 x 24 241 x 251 x 22	JG251R	406 x 253 x 24 404 x 251 x 22	JG211R	485 x 253 x 24 483 x 251 x 22	JG221R	540 x 253 x 24 538 x 251 x 22
 44	JG112R	243 x 253 x 44 241 x 251 x 42	JG252R	406 x 253 x 44 404 x 251 x 42	JG212R	485 x 253 x 44 483 x 251 x 42	JG222R	540 x 253 x 44 538 x 251 x 42
 64	JG113R	243 x 253 x 64 241 x 251 x 62	JG253R	406 x 253 x 64 404 x 251 x 62	JG213R	485 x 253 x 64 483 x 251 x 62	JG223R	540 x 253 x 64 538 x 251 x 62
 94	JG114R	243 x 253 x 94 241 x 251 x 92	JG254R	406 x 253 x 94 404 x 251 x 92	JG214R	485 x 253 x 94 483 x 251 x 92	JG224R	540 x 253 x 94 538 x 251 x 92

Basket Inserts



		Outside dim. Inside dim.		Outside dim. Inside dim.		Outside dim. Inside dim.
 34	JF110R	196 x 252 x 34 194 x 237 x 33	JF250R	305 x 252 x 34 303 x 237 x 33	JF210R	438 x 252 x 34 436 x 237 x 33

Lid for Baskets JF, JG



	 Standard 1/2	Outside dim.	 Standard 3/4	Outside dim.	 Standard 1/1 short	Outside dim. Inside dim.	 Standard 1/1 long	Outside dim. Inside dim.
 JF117R		247 x 257 x 18	JF257R		410 x 257 x 18	JF217R		489 x 257 x 18
						JF227R		544 x 257 x 18
						JH217R		489 x 257 x 18 with 8 mm perforation

Aesculap JF/JG baskets are made of stainless steel 1.4301 and of PEEK (feet, closure system basket lid)

All dimensions: Length x Width x Height (mm)

Baskets

Special trays



4 mm perforation	Outside dim. Inside dim.
JF155R (lid included)	236 x 72 x 27 230 x 66 x 20



Lid, fitting JF156R and JF155R

4 mm perforation	Outside dim. Inside dim.
JF157R	236 x 72 x 11




For small objects, two trays fitting into an AESCULAP® mini sterile container, 2 each fitting into JK188 or JN188

Mini	Outside dim. Inside dim.
JF149R (2 mm perforation)	135 x 120 x 47 128 x 113 x 38



Tray base

4 mm perforation	Outside dim. Inside dim.
JF156R	232 x 68 x 61 227 x 61 x 54

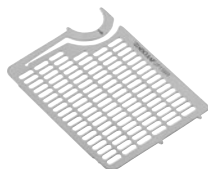

Tray divider, fitting JF156R

JF158R	



For small objects or micro instruments

4 mm perforation	Outside dim. Inside dim.
JF150R	235 x 137 x 42 229 x 131 x 35



For micro instruments fitting into an AESCULAP® mini sterile container, fitting into JK188 or JN188

Mini	Outside dim. Inside dim.
JF159R (4 mm perforation)	267 x 125 x 49 261 x 118 x 41



suitable sterile containers see page 43



suitable sterile containers see page 44

Special trays



Tray with lid

Dental	Outside dim.
	Inside dim.
JG487R	188 x 84 x 30 181 x 82 x 23



Tray with lid

Dental	Outside dim.
	Inside dim.
JG488R	188 x 125 x 30 181 x 123 x 23



For small objects

	Outside dim.
	Inside dim.
JF145R 2 mm perforation	85 x 44 x 25 80 x 39 x 18
JF146R 4 mm perforation	85 x 44 x 25 80 x 39 x 18



For small objects

4 mm perforation	Outside dim.
	Inside dim.
JF140R	80 x 80 x 33 75 x 75 x 26



For small objects

4 mm perforation	Outside dim.
	Inside dim.
JF147R	354 x 44 x 35 349 x 39 x 28



For small objects

4 mm perforation	Outside dim.
	Inside dim.
JF148R	145 x 44 x 25 140 x 39 x 18

Baskets

Special trays



Tray without inserts, with lid handle fitting into dental sterile container JN294

	Outside dim. Inside dim.
JF264R	273 x 176 x 41 267 x 170 x 34



Tray with lid fitting into dental sterile container JN295

	Outside dim. Inside dim.
JG285R	273 x 176 x 60 267 x 170 x 53



Tray without lid fitting into dental sterile container JN296

	Outside dim. Inside dim.
JF282R	269 x 171 x 93 267 x 169 x 91



Lid for tray JF282R

	Outside dim.
JF289R	273 x 176 x 16



Tray, contains three instrument racks (JG392), with lid handle, fitting into dental sterile container JN294

	Outside dim. Inside dim.
JG380R	273 x 84 x 41 267 x 78 x 34



Tray, contains three instrument racks (JG391), with lid handle, fitting into dental sterile container JN294

	Outside dim. Inside dim.
JG381R	273 x 176 x 30 267 x 170 x 23

Special trays



Tray, with perforation for retractor, contains three instrument racks (JG391), with lid handle, fitting into dental sterile container JN294

	Outside dim.
	Inside dim.
JG383R	273 x 176 x 30 267 x 170 x 23



Tray, contains three instrument racks (JG393), with lid handle, fitting into dental sterile container JN294

	Outside dim.
	Inside dim.
JG384R	273 x 176 x 41 267 x 170 x 34



Tray, without inserts, with lid handle, fitting into dental sterile container JN294

	Outside dim.
	Inside dim.
JG386R	273 x 176 x 41 267 x 170 x 34



Tray, contains the instrument racks JG400 and JG401 for 6 instruments

	Outside dim.
	Inside dim.
JG387R	188 x 84 x 30 181 x 78 x 23



Tray, contains the instrument racks JG400 and JG401 for 9 instruments


	Outside dim.
	Inside dim.
JG388R	188 x 125 x 30 181 x 119 x 23



Tray, contains five instrument racks (2 x JG394, 2 x JG400, 1 x JG401), with lid handle, fitting into dental sterile container JN294


	Outside dim.
	Inside dim.
JG389R	273 x 176 x 30 267 x 170 x 23

Scope trays

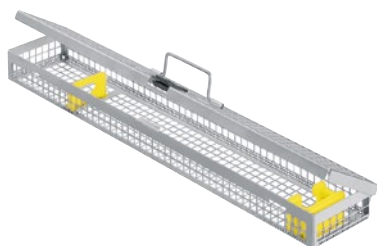


Tray for one scope of sheath length up to 340 mm.

	Outside dim.
JF431R	454 x 84 x 41 448 x 78 x 34

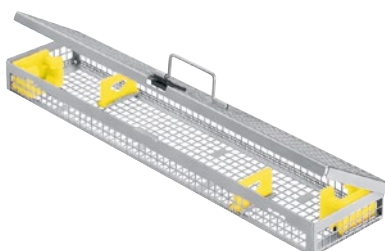


Tray for two scopes of sheath length up to 340 mm.

	Outside dim. Inside dim.
JF432R	454 x 104 x 41 448 x 98 x 34

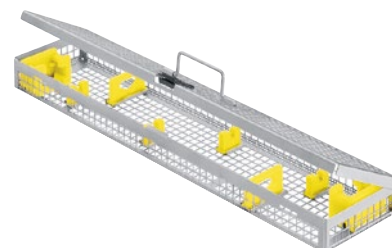

Tray for two scopes of sheath length up to 340 mm, with additional support.

	Outside dim. Inside dim.
JF433R	454 x 104 x 41 448 x 98 x 34

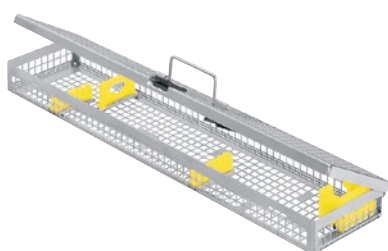

Tray for mini hysteroscope of < 4 mm diameter


	Outside dim. Inside dim.
JF434R	454 x 104 x 41 448 x 98 x 34



Tray for one arthroscope of 2.7 mm or 4 mm sheath diameter and sheath length 140 - 200 mm.


	Outside dim. Inside dim.
JF436R	273 x 84 x 41 267 x 78 x 34

Scope trays

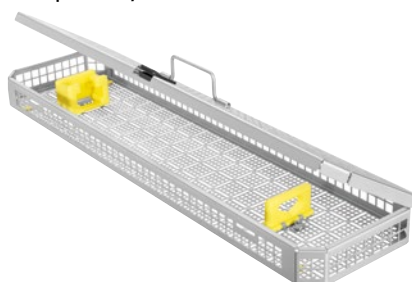

Tray like JF435R prepared for one scope and one light conducting cable (incl. 1x JG367 + 1x JG368), fitting scope container see page 45 ff.

	Outside dim. Inside dim.
JF428R	500 x 132 x 44 494 x 124 x 30

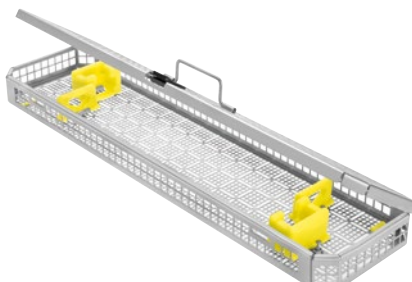

Tray like JF435R prepared for two scopes and one light conducting cable (incl. 2x JG367 + 2x JG368), fitting scope container see page 45 ff.


	Outside dim. Inside dim.
JF429R	500 x 132 x 44 494 x 124 x 30



Tray for up to two laparoscopic scopes and one light conducting cable (scope holders JG367 + JG368 please order separately), fitting scope container see page 45 ff.


	Outside dim. Inside dim.
JF435R	500 x 132 x 44 494 x 124 x 30

Storage rack for MIC instruments and trocars

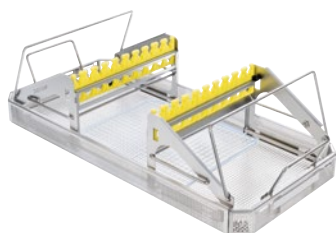

Storage rack set for MIS instruments and trocars, consisting of: Storage rack PL963R, instrument basket JF522R and silicone mat JF938, fitting JK444 and JK489.

Standard 1/1	Outside dim. Inside dim.
PL960R	540 x 253 x 166

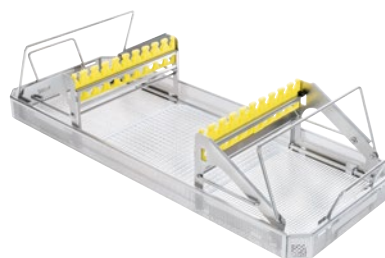

Storage rack set for MIS instruments, XXL, consisting of: Storage rack PL964R, instrument basket JF532R and silicone mat JF938, for instruments with shaft length up to 420 mm, fitting JN445 and JK490.

XXL	Outside dim. Inside dim.
PL961R	660 x 253 x 166

Storage racks



Storage rack, standard


	Outside dim.
PL963R	487 x 250 x 153



Storage rack, XXL


	Outside dim. Inside dim.
PL964R	587 x 250 x 153

Storage frame

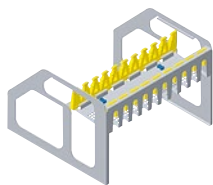

Storage frame
for bone punches


	Outside dim. Inside dim.
JF120R	240 x 225 x 151

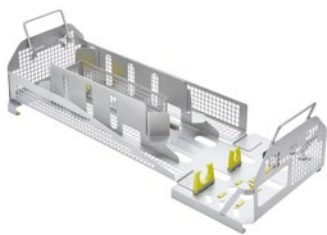

Da Vinci SP* Camera rack to be used
with the Endocontainer JN445 bottom
and JK490 lid

	Outside dim. Inside dim.
JF449R	660 x 253 x 163

*Registered trademark by Intuitive Surgical

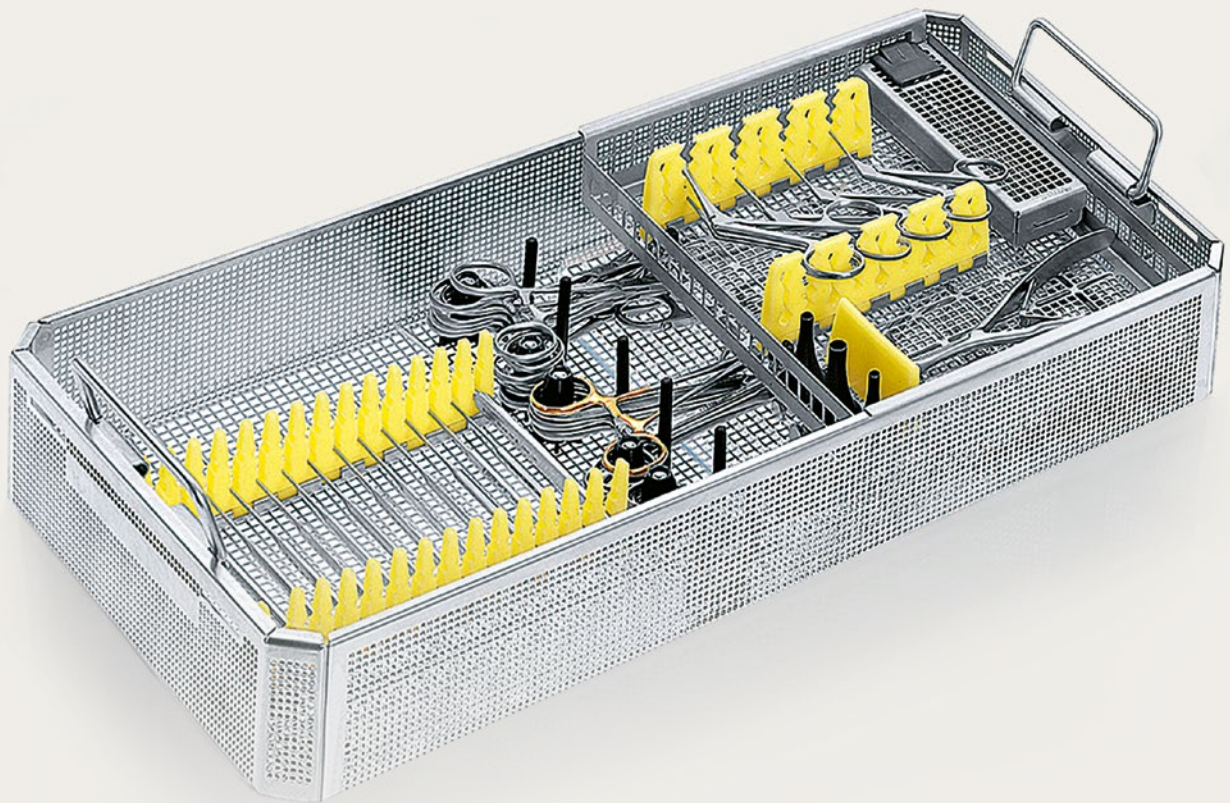

Aesculap InOS® – Instrument Organization System

Aesculap InOS[®] – Instrument Organization System


Organization and protection of surgical instruments

Surgical instruments are valuable assets for every hospital, that is why they should be processed as efficiently and gently as possible.


Therefore, protecting the instruments during sterilisation, as well as in transport to sterile supply storage and to the operating room, is of crucial importance.

The Aesculap organization system for instrument trays facilitates the storage especially of delicate instruments and small items.


Transsphenoidal
hypophysectomy Set

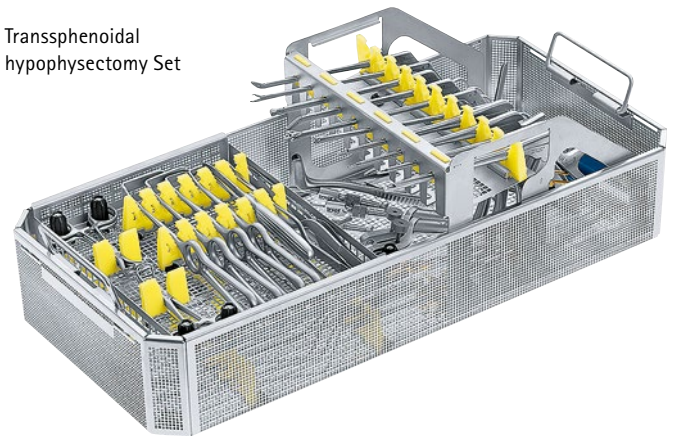


Tympanoplasty Set

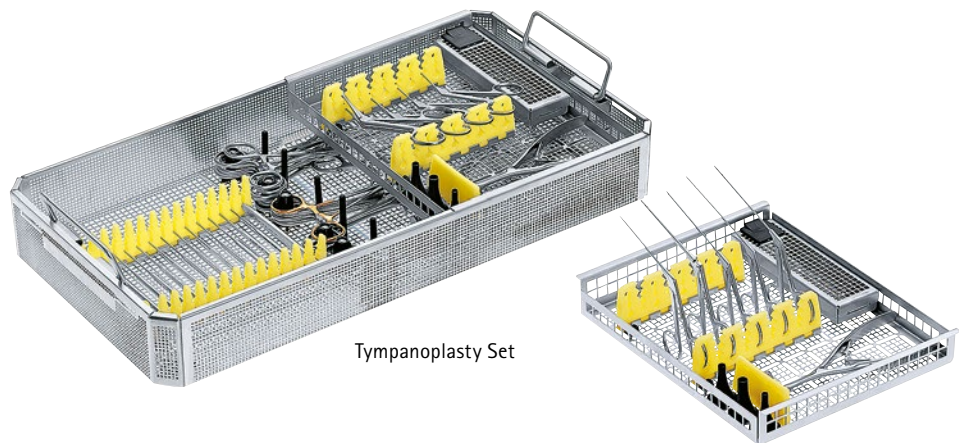


Skull Set

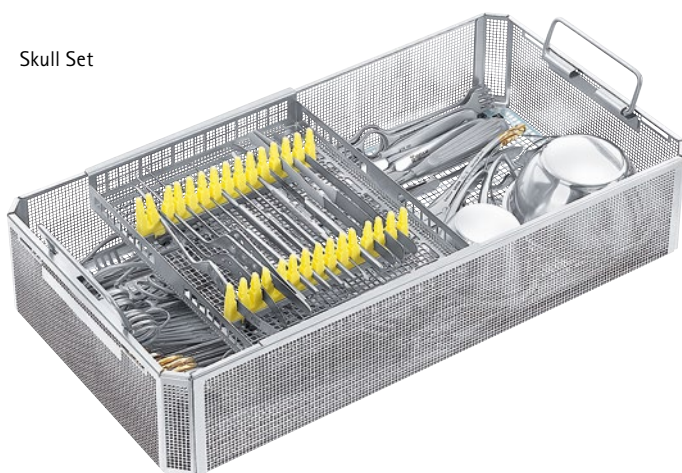

Laminectomy Set

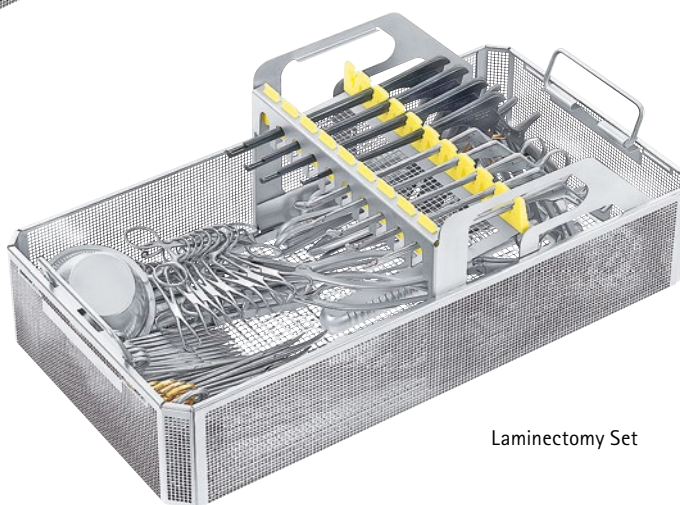


CABG Micro accessory
(Coronary Artery Bypass Graft)

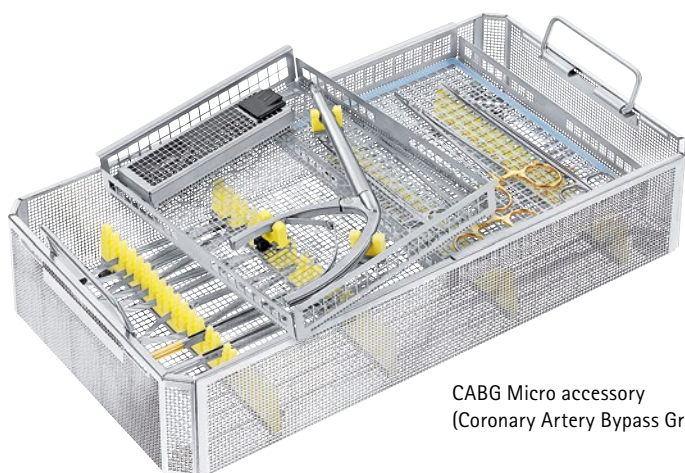


Aesculap InOS® – Instrument Organization System

Instrument racks and holders

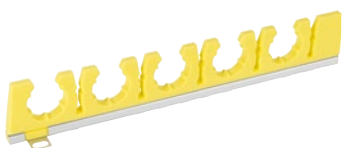Instrument racks* (mounting A)			Interspaces	Outside dim.	Art. no.	
			7	120 x 40 mm	JG304	
			9	160 x 40 mm	JG305	
			15	240 x 40 mm	JG306	
			9	120 x 40 mm	JG308	
			11	160 x 40 mm	JG309	
			17	240 x 40 mm	JG310	
			6	120 x 40 mm	JG312	
			8	160 x 40 mm	JG313	
			12	240 x 40 mm	JG314	
			10	120 x 40 mm	JG316	
			14	160 x 40 mm	JG317	
			20	240 x 40 mm	JG318	
			20	120 x 40 mm	JG327	
			for curettes	10	240 x 41 mm	JG311
Instrument rack for bipolar forceps or other angled instruments			Interspaces	Outside dim.	Art. no.	
			29	234 x 41 mm	JG319	


Instrument racks for micro instruments* (mounting A)			Interspaces	Outside dim.	Art. no.
			20	120 x 21 mm	JG320
			7	123 x 18 mm	JG321
			6	123 x 18 mm	JG322
			13	123 x 18 mm	JG323
			7	123 x 45 mm	JG324

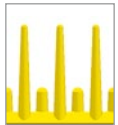Instrument racks for micro instruments* (mounting A)			Interspaces	Outside dim.	Art. no.
		Positioning rack	19	160 x 30 mm	JG325
			29	240 x 30 mm	JG326
		Support for spring type micro instruments – handle end	4	120 x 40 mm	JG328
			2	60 x 40 mm	JG333
		Support for spring type micro instruments – jaw end	4	120 x 40 mm	JG329
			2	60 x 40 mm	JG334

Aesculap InOS® – Instrument Organization System

Holders with rail* (mounting A)	Diameter	Art. no.
	5 mm	JG340
	7.5 mm	JG341
	8.5 mm	JG342
	10 mm	JG343
	12 mm	JG344

Metal clamp with silicone inlay (mounting A)	Diameter	Art. no.
	15 - 20 mm	JG352R
	20 - 30 mm	JG353R
	30 - 45 mm	JG354R

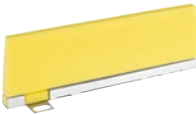
Separating racks* (mounting A)	Length	Art. no.
	120 x 40 mm	JG358
	160 x 40 mm	JG359
	243 x 40 mm	JG360


Note | Unless noted otherwise, please order the screws JG303R separately.

Aesculap InOS® – Instrument Organization System

Holders (mounting A)		Diameter	Art. no.	Diameter	Art. no.
		15 mm	JG330	10 mm	JG335
		25 mm	JG331	7.5 mm	JG336
		30 mm	JG332	5 mm	JG337

Scope holders* (mounting A)		Interspaces	Art. no.
		Eyepiece holder	JG367
		Sleeve holder	JG368
		Support for angled optics	JG369

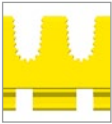Racks for basket inlays (mounting C)	Interspaces	Outside dim.	Art. no.
 		240 x 30 mm	JG370
 	10	234 x 30 mm	JG371
 	10	227 x 30 mm	JG372

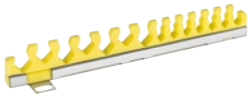Racks for special trays (mounting A)	Interspaces	Outside dim.	Art. no.
	6	76 x 33 mm	JG396
	14	160 x 23 mm	JG397 (incl. holding screws)
	13	160 x 34 mm	JG398 (incl. holding screws)

Aesculap InOS® – Instrument Organization System


Instrument racks for special trays with 30 mm height (mounting C)	Interspaces	Outside dim.	Art. no.
	13	159 x 14 mm	JG400
	6	74 x 22 mm	JG394
	10	115 x 22 mm	JG395
	14	158 x 22 mm	JG391

Note | Unless noted otherwise, please order the screws JG303R separately.

Instrument racks for special trays with 41 mm height (mounting C)	Interspaces	Outside dim.	Art. no.
	6	76 x 33 mm	JG392
	13	158 x 33 mm	JG393
		158 x 30 mm	JG390

Top holder		Outside dim.	Art. no.
		113 x 21 mm	JG378
		113 x 31 mm	JG379
	Press pad, for suppressing the instrument, to assemble to the lower side of the lid.	150 x 15 mm	JG399
		158 x 14 mm	JG401


Aesculap InOS® – Instrument Organization System




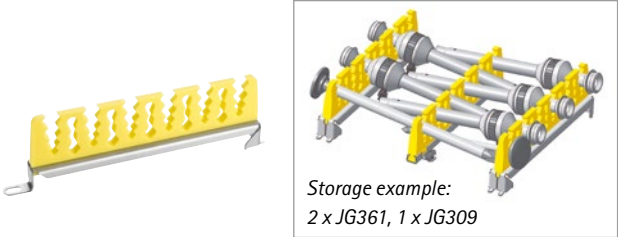
Trocar rack (mounting A)	Interspaces	Outside dim.	Art. no.
 Storage example: 2 x JG361, 1 x JG309	11	167 x 53 mm	JG361


Fixation elements for mounting A		Art. no.
	Screw set, to be used with GB459R. Order unit (PAK): Pack of 10	JG303R
	Eccos & Aesculap InOS® mounting kit	GB459R

Fixation pins for ping handle instruments	Diameter	Art. no.
	25 mm, height 29 mm	JF346
	6 mm, height 47 mm	JF343*

***Note** | does not fit into JJ baskets
(AESCULAP Aicon®)


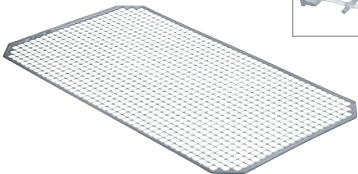
Silicone mats (Mesh)		Outside dim.	Art. no.
 	fitting to 1/2 basket	240 x 250 mm	JF938
	fitting to 3/4 basket	402 x 250 mm	JF940
	fitting to 1/1, short basket	480 x 250 mm	JF939
	fitting to 1/1, long basket	536 x 250 mm	JF941

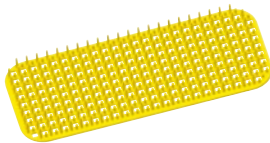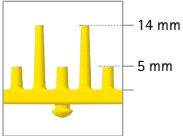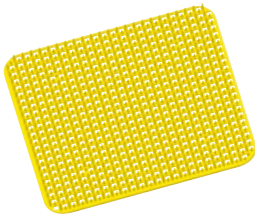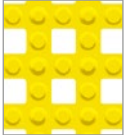Silicone mats with perforation (Nubs)		Outside dim.	Art. no.
 	fitting to JF159R	248 x 102 mm	JF944
	fitting to 1/2 basket	248 x 237 mm	JF945
	fitting to 3/4 basket	394 x 242 mm	JF947
	fitting to 1/1, short basket	472 x 242 mm	JF948
 	fitting to 1/1, long basket	517 x 242 mm	JF949

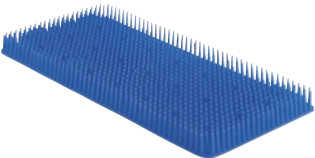Positioning mat		Outside dim.	Art. no.
		470 x 230 x 30 mm	JF932
 		276 x 125 x 18 mm	JF934

Aesculap InOS® – Instrument Organization System

Storage solution for Hegar uterine dilators

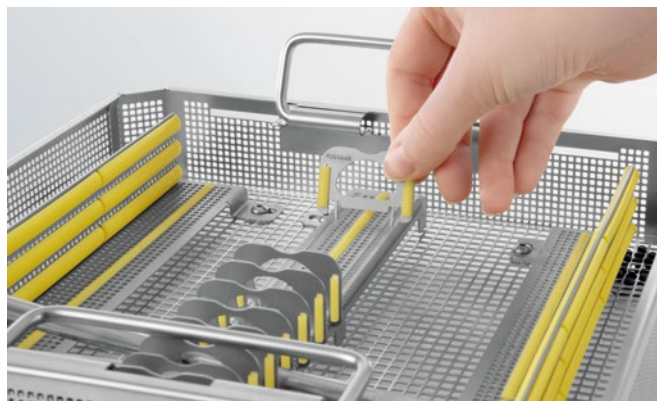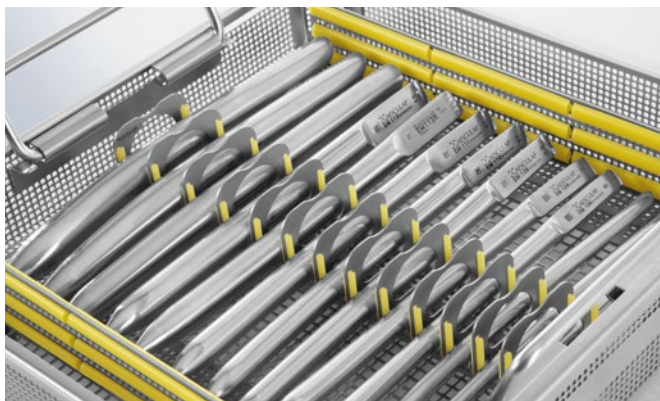


Example:

1 x JG362R center bar
2 x JG363R angle bar
10 to 15 dividers JG364R
1 x JF112R 1/2 basket
1 x JF117R 1/2 basket lid (not depicted)




Also fitting

1 x JJ262R AESCULAP Aicon® 1/2 basket without lid or
1 x JJ263R AESCULAP Aicon® 1/2 basket with lid
JJ267R (not depicted)

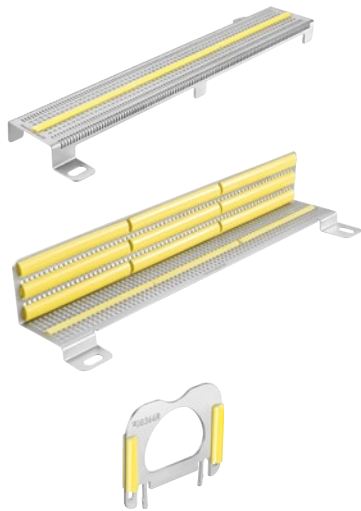Storage solution for uterine dilators		Outside dim.	Art. no.
	Center bar with fixation screws	213 x 67 x 16 mm	JG362R
	Angle bar with fixation screws	213 x 50 x 42 mm	JG363R
	Divider for center bar JG362R	40 x 5 x 30 mm	JG364R

Mounting instructions

Mounting A for instrument racks to be fixed with screws

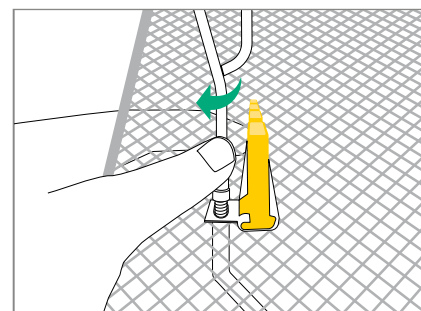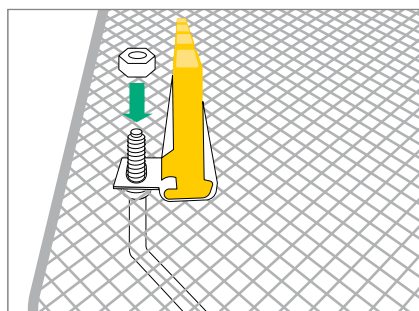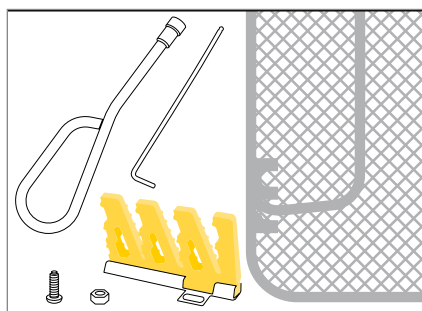




Materials required:

- Screw set JG303R
- Eccos Et Aesculap InOS® mounting kit GB459R
- Instrument holder such as JG341
- AESCULAP baskets (e. g. AESCULAP Aicon® JJ series or JF/JG series)

Place the instrument holder in the required position in the AESCULAP basket, making sure that the holder's fixation lug fit onto the basket's lattice-work correctly. Push the screw through from below and fix it with the nut from above.

Hold the screw with the hex-wrench, put on the socket wrench with the other hand and tighten the nut-bolt connection firmly. Make sure that the bolt head only projects slightly above the basket bottom after installation.

Mounting C for silicone racks in basket inserts

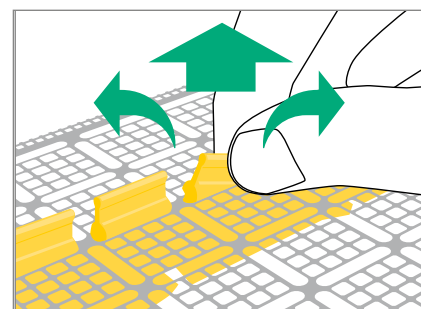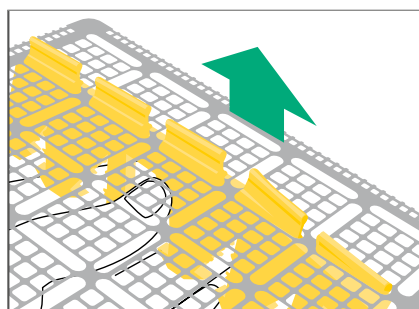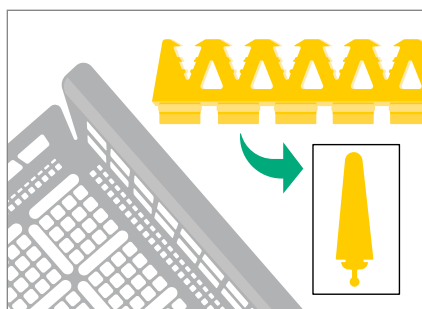

Materials required:

- Basket inserts JF110R, JF250R or JF210R
- Instrument rack, such as JG372

Push the fixation lugs through the longitudinal slots in the basket.

To pull the lugs through the slot completely, move left and right while applying tension.

Aesculap InOS® – Instrument Organization System

Holding clamps

Not for AESCULAP Aicon® baskets
(JJ-series)



	Outside dim. Effective dim.	
JG691R	24 x 28 x 25 Ø 7-12 mm	Holding clamp



	Outside dim. Effective dim.	
JG692R	24 x 28 x 29 Ø 12-16 mm	Holding clamp



	Outside dim. Effective dim.	
JG693R	24 x 28 x 31 Ø 16-26 mm	Holding clamp



	Outside dim. Effective dim.	
JG694R	24 x 28 x 40 Ø 16-36 mm	Holding clamp



	Outside dim. Effective dim.	
JG695R	24 x 28 x 44 Ø 36-46 mm	Holding clamp

Fixation elements

Not for AESCULAP Aicon® baskets
(JJ-series)



	External dim.
JG696R	Ø 6.0 x 20
JG702R	Ø 4.9 x 20
JG703R	Ø 4.5 x 22
JG706R	Ø 6.0 x 28
JG708R	Ø 6.0 x 45

Holding pin – e.g. for motor accessories

Holding pin – e.g. for bone flaps

Holding pin – e.g. for coriotome heads

Holding pin – e.g. for motor accessories

Holding pin – e.g. for motor accessories

Aesculap InOS® – Instrument Organization System

Stop frames



	Outside dim.
JF921R	110 x 15 x 28

Stop frame,
with fastening screws.



	Outside dim.
JF922R	174 x 15 x 28

Stop frame,
with fastening screws.



	Outside dim.
JF923R	244 x 15 x 28

Stop frame,
with fastening screws.



	Outside dim.
JF378R	164 x 27 x 12

Stop frame,
with fastening screws.



Accessories

Accessories

Identification labels



JG785B	red
JG786B	blue
JG787B	green
JG788B	gold
JG790B	black
JG792B	olive-brown
JG793B	violet
JG794B	orange

Identification labels with lettering one line or two lines (character height 5 mm), up to 13 characters per line. Dimensions: 56 x 18 mm for all Standard sterile containers except JK299, JK309 and JK339.

Note: Process aluminum identification labels with pH-neutral cleaning agent.



JG785U	red
JG786U	blue
JG787U	green
JG788U	gold
JG790U	black
JG792U	olive-brown
JG793U	violet
JG794U	orange

Identification labels without lettering

Specifications as above

Note: Process aluminum identification labels with pH-neutral cleaning agent.



JG645B	red
JG646B	blue
JG647B	green
JG648B	gold
JG650B	black
JG652B	olive-brown
JG653B	violet
JG654B	orange

Identification labels with lettering one line or two lines (character height 5 mm), up to 13 characters per line. Dimensions: 50 x 18 mm for mini-, dental-, scope sterile container, disposal container and JK299, JK309, JK339.

Note: Process aluminum identification labels with pH-neutral cleaning agent.



JG645U	red
JG646U	blue
JG647U	green
JG648U	gold
JG650U	black
JG652U	olive-brown
JG653U	violet
JG654U	orange

Identification labels without lettering

Specifications as above

Note: Process aluminum identification labels with pH-neutral cleaning agent.



Identification labels



	Outside dim.
JG736	56 x 18

Bar code labels, plastic, general, suitable for the label holder on AESCULAP® sterile containers. type 128, 10 characters max. code length, 10 characters max. lettering length.

Special note:
Please specify the required designation in your order. Bar code and designation represent the same content. Also fitting into holder JG730.



	Outside dim.
IC-V006	56 x 18

Bar code labels, plastic, for instacount®, fitting into the label holder on AESCULAP® sterile containers. type 128, 10 characters max. code length, 20 characters max. lettering length.

Special note:
Please specify the required designation and the serial number of the instacount set.



	Outside dim.
JG730U	50 x 35 x 13

Holder for identification labels, for mounting on the trays, blank.



	Outside dim.
JG730B	50 x 35 x 13

Holder for identification labels, for mounting on the trays. Lettered (one line, max. 14 characters/line, character height max. 4 mm).

Accessories

Identification labels



JG625R

Outside dim.

50 x 20 x 4

Holder for identification labels,
to be mounted on baskets.



JG636R

Outside dim.

33 x 20

Identification label for baskets, with lettering
(one-line or two-line), up to 8 characters
per line, character height 5 mm.



JG616U

Outside dim.

58 x 23

Identification label for baskets, for barcode tags.

JG616B

58 x 23

Identification label for baskets, with lettering
(one-line or two-line), up to 14 characters per
line, character height 5 mm.



JG617U

Outside dim.

89 x 30

Identification label for baskets, for barcode tags.

JG617B

89 x 30

Identification label for baskets, with lettering
(one-line or two-line), up to 20 characters per
line, character height 5 mm.



JJ586U

Outside dim.

96 x 34

Identification label for baskets, for barcode tags.

JJ586B

96 x 34

Identification label for baskets, with lettering
(one-line, two-line or three-line), up to 13
characters per line, character height 5 mm.



JJ587U

Outside dim.

105 x 54

Identification label for baskets, for barcode tags.

JJ587B

105 x 54

Identification label for baskets, with lettering
(one-line, two-line or three-line), up to 13
characters per line, character height 5 mm.

Disposal stamps



	Outside dim.
JG673B	40 x 35 x 9
JG673U	40 x 35 x 9

Disposal stamp with lettering
(one-line or two-line), up to 14 characters per line.

Disposal stamp, blank
Order unit: ST = pieces



	Outside dim.
JG675B	40 x 35 x 9
JG675U	40 x 35 x 9

Disposal stamp with lettering
one line, up to 3 digits (numerical only).

Disposal stamp, blank
Order unit: PR = pair



	Outside dim.
JG690R	393 x 85 x 155

Storage stand
for disposal stamps.

Fixation clamps



	Outside dim.
JF484R	65 x 20 x 20

Fixation clamp
for clamping together baskets and containers
that are arranged next to each other on the prep
table.



	Outside dim.
JF509R	50 x 20 x 12

Fixation clamp
for wrapping drapes.

Accessories

Filter accessories



JK100

Universal filter retention plate for reusable filters and paper filters, fitting for all lid and body components of standard sterile containers.



JG731R

Bar code identification ring for filter retention plate. For recording the cycle count of reusable filter JK090. Lettering acc. to specification, up to 6 digits with bar code BC128.



JK098

Outside dim.

245 x 100 x 6

Universal filter retention plate for reusable and paper filters, fitting for all lid and base components of mini, dental and scope sterile containers.



JK096

Outside dim.

197 x 197 x 40

Filter dispenser for round filters.



JP050

Filter insert for PRIMELine® and PRIMELine® PRO lids validated for 5,000 sterilization cycles.



JP055

Mounting handle for filter insert, PRIMELine® and PRIMELine® PRO lids.

Consumables



Consumables

Reusable filters



Standard XXL	Outside dim.
JK090	Ø 190

AESCULAP® reusable filter, round for use in up to 1,000 cleaning and sterilization cycles.

Order unit (PAK): Pack of 10



Standard	Outside dim.
JK089	225 x 225

AESCULAP® reusable filter, square for legacy standard sterile containers, service life 1,000 cleaning and sterilization cycles.

Order unit (PAK): Pack of 10



Mini Dental Scope	Outside dim.
JK091	95 x 215

AESCULAP® reusable filter for mini, dental and scope sterile containers with filter retention plate JK098, service life 1,000 cleaning and sterilization cycles.

Order unit (PAK): Pack of 2



Dental	Outside dim.
JK092	171 x 231

AESCULAP® reusable filter for dental sterile containers without filter retention plate JK098, service life 1,000 cleaning and sterilization cycles.

Order unit (PAK): Pack of 2

Paper filters



Standard	Outside dim.
XXL	
JK095	Ø 190

AESCULAP® paper filter
with process indicator for steam sterilization*.

Order unit (PAK): Pack of 100



Mini	Outside dim.
Dental	
Scope	
JF166	118 x 235

AESCULAP® paper filter for mini, dental and
scope sterile containers with process
indicator for steam sterilization*, for use
with filter retention plate JK098.

Order unit (PAK): Pack of 100



Dental	Outside dim.
JF066	171 x 232

AESCULAP® paper filter for dental sterile
containers with process indicator for steam
sterilization*, for use without filter retention
plate JK098.


Order unit (PAK): Pack of 100



Standard	Outside dim.
JF481	230 x 230

AESCULAP® paper filter, square
for legacy standard sterile containers without
process indicator.

Order unit (PAK): Pack of 100

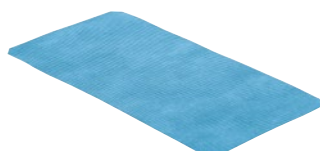

Scope	Outside dim.
JF167	235 x 118

Disposable filter for the hydrogen peroxide
sterilisation in Sterrad® 100S, Sterrad® 200,
Sterrad® NX, Sterrad® 100NX, Steris V-Pro® 1
and Steris V-Pro® 1 Plus.





Order unit (PAK): Pack of 100

* acc. to EN ISO 11140-1 type 1

All dimensions: Length x Width x Height (mm)

Consumables

Indicator labels

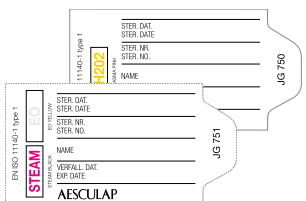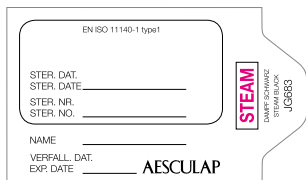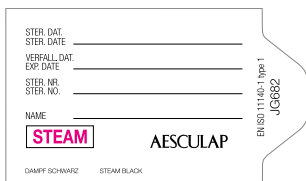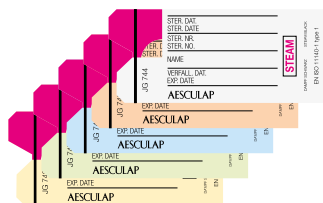

Standard	Color	Outside dim.
JG744	white	80 x 35
JG745	red	80 x 35
JG746	blue	80 x 35
JG747	green	80 x 35
JG748	yellow	80 x 35

Standard	Outside dim.
JG743	80 x 35

Mini Dental	Outside dim.
JG687	60 x 18

	Outside dim.
JG682	67 x 40

	Outside dim.
JG683	67 x 40

Scope	Outside dim.
JG750	67 x 35
JG751	67 x 35

Indicator seal with adhesive back and process indicator for steam sterilization*; safeguarding the sterile container against unauthorized opening. For standard sterile containers.
Order unit (PAK): Pack of 1,000

Seal without adhesive back and without process indicator for steam sterilization; safeguarding the sterile container against unauthorized opening. For standard sterile containers.
Order unit (PAK): Pack of 1,000

Indicator label with adhesive back and process indicator for steam sterilization* for mini and dental sterile containers, container bottom JK339.
Order unit (PAK): Pack of 100

Indicator label with process indicator for steam sterilization* for legacy standard sterile containers.
Order unit (PAK): Pack of 1,000

Indicator label with adhesive back and process indicator for steam sterilization, for legacy standard sterile containers.
Order unit (PAK): Pack of 1,000

Identification label with adhesive back and process indicator for sterilisation with vaporized hydrogen peroxide*.
Order unit (PAK): Pack of 100

Identification label with adhesive back and process indicator for steam sterilisation* and sterilisation with ethylene oxide*.
Order unit (PAK): Pack of 100

Seals

	Mini Dental Standard	JG738 JG739	Seal - safeguarding the sterile container against unauthorized opening. Order unit (PAK): Pack of 100 Order unit (PAK): Pack of 1,000
	Mini Dental Standard	JG740	Repair stamp - for marking instruments for repair or sterile containers holding materials due for repair. Order unit (PAK): Pack of 1,000
	Mini Dental Standard	JG752 JG753	Seal - safeguarding the sterile container against unauthorized opening, with process indicator for steam sterilization*. Order unit (PAK): Pack of 100 Order unit (PAK): Pack of 1,000

* acc. to EN ISO 11140-1 type 1

Consumables

Maintenance products



JG598

STERILIT® I
AESCULAP® drip oiler for lubricating the hinges of the container locking mechanism, instrument joints etc. prior to sterilization.

Contents: 50 ml



JG600

STERILIT® I – including 6 spraying nozzles
AESCULAP® oil spray for lubricating the hinges of the container locking mechanism, instrument joints etc. prior to sterilization.

Contents: 300 ml
Order unit: Pack of 6



JG601

Eloxal cleaner
for cleaning and care of sterile container aluminium surfaces.

Contents: 300 g

Labels, bar code labels



IC-V004

Roll of labels, double self-adhesive, outward rolled, without indicator, for cutting tool.

Length 70 m (approx. 2,000 labels 35 mm high), width 75 mm, roll core diameter 3", order unit (ST).



IC-V001

Roll of labels, double self-adhesive, outward roll, with process indicator for steam sterilization acc. to EN ISO 11140-1 type 1, for cutting instruments.

Length 70 m (approx. 2,000 labels 35 mm high), width 77 mm, roll core diameter 3", order unit (ST).



IC-V002

Individual labels, double self-adhesive, outward rolled, without indicator, for dispenser sensor.

Label height 35 mm, label width 75 mm, roll core diameter 3", order unit (ST).



IC-V005

Individual labels, double self-adhesive, outward rolled, for dispenser sensor.

Label height 35 mm, label width 75 mm, roll core diameter 3", order unit (ST).



IC-V003

Thermotransfer ribbon, wax, inside print.

Length 360 m, width 76.5 mm, core diameter 1", order unit (ST).



IC-V008

Thermotransfer ribbon, wax, outside print.

Length 300 m, width 80 mm, core diameter 1", order unit (ST).

Consumables

Labels, printer accessories



IC-V012

Individual labels, double side self adhesive, outside wound, with process indicator, with EN ISO 11140-1 type 1, for cutting tool.

1.500 labels perforated, label height 35 mm, width 78 mm; core diameter 3".



IC-V013

Endless labels, single side self adhesive, outside wound, with process indicator for steam sterilization in compliance with EN ISO 11140-1 type 1, for cutting tool.

Length 70 m (approx. 2000 labels at a height of 35 mm); width 77 mm; core diameter 3".



IC-V017

Tag label, special foil, triple-layer, white.

Outer format 75 x 265 mm
Order unit (PAK): Pack of 500

Index

Index

Article No.	Page	Article No.	Page	Article No.	Page
GB459R	74	JF221R	54	JF944	75
IC-V001	93	JF222R	54	JF945	75
IC-V002	93	JF223R	54	JF947	75
IC-V003	93	JF224R	54	JF948	75
IC-V004	93	JF227R	55	JF949	75
IC-V005	93	JF232R	44		
IC-V006	83	JF236R	44	JG111R	54
IC-V008	93	JF250R	55	JG112R	54
IC-V012	94	JF251R	54	JG113R	54
IC-V013	94	JF252R	54	JG114R	54
IC-V017	94	JF253R	54	JG211R	54
		JF254R	54	JG212R	54
JF066	89	JF257R	55	JG213R	54
JF110R	55	JF264R	58	JG214R	54
JF111R	54	JF282R	58	JG221R	54
JF112R	54	JF289R	58	JG222R	54
JF113R	54	JF343	74	JG223R	54
JF114R	54	JF346	74	JG224R	54
JF117R	55	JF378R	80	JG251R	54
JF120R	62	JF428R	61	JG252R	54
JF140R	57	JF429R	61	JG253R	54
JF145R	57	JF431R	60	JG254R	54
JF146R	57	JF432R	60	JG285R	58
JF147R	57	JF433R	60	JG303R	74
JF148R	57	JF434R	60	JG304	66
JF149R	56	JF435R	42, 61	JG305	66
JF150R	56	JF436R	60	JG306	66
JF155R	56	JF449R	45, 62	JG308	66
JF156R	56	JF481	89	JG309	66
JF157R	56	JF484R	85	JG310	66
JF158R	56	JF509R	85	JG311	66
JF159R	56	JF921R	80	JG312	66
JF166	42, 89	JF922R	80	JG313	66
JF167	43, 89	JF923R	80	JG314	66
JF210R	55	JF932	75	JG316	66
JF211R	54	JF933	44	JG317	66
JF212R	54	JF934	75	JG318	66
JF213R	54	JF938	75	JG319	66
JF214R	54	JF939	75	JG320	67
JF217R	55	JF940	75	JG321	67
		JF941	75	JG322	67

Article No.	Page	Article No.	Page	Article No.	Page
JG323	67	JG384R	59	JG653U	82
JG324	67	JG386R	59	JG654B	82
JG325	67	JG387R	59	JG654U	82
JG326	67	JG388R	59	JG673B	85
JG327	66	JG389R	59	JG673U	85
JG328	67	JG390	73	JG675B	85
JG329	67	JG391	72	JG675U	85
JG330	70	JG392	73	JG682	90
JG331	70	JG393	73	JG683	90
JG332	70	JG394	72	JG687	90
JG333	67	JG395	72	JG690R	85
JG334	67	JG396	71	JG691R	78
JG335	70	JG397	71	JG692R	78
JG336	70	JG398	71	JG693R	78
JG337	70	JG399	73	JG694R	78
JG340	68	JG400	72	JG695R	78
JG341	68	JG401	73	JG696R	79
JG342	68	JG487R	57	JG702R	79
JG343	68	JG488R	57	JG703R	79
JG344	68	JG598	92	JG706R	79
JG352R	69	JG600	92	JG708R	79
JG353R	69	JG601	92	JG730B	83
JG354R	69	JG616B	84	JG730U	83
JG358	69	JG616U	84	JG731R	86
JG359	69	JG617B	84	JG736	83
JG360	69	JG617U	84	JG738	91
JG361	74	JG625R	84	JG739	91
JG362R	76	JG636R	84	JG740	91
JG363R	76	JG645B	82	JG743	90
JG364R	76	JG645U	82	JG744	90
JG367	70	JG646B	82	JG745	90
JG368	70	JG646U	82	JG746	90
JG369	70	JG647B	82	JG747	90
JG370	71	JG647U	82	JG748	90
JG371	71	JG648B	82	JG750	43, 90
JG372	71	JG648U	82	JG751	42, 90
JG378	73	JG650B	82	JG752	91
JG379	73	JG650U	82	JG753	91
JG380R	58	JG652B	82	JG785B	82
JG381R	58	JG652U	82	JG785U	82
JG383R	59	JG653B	82	JG786B	82

Index

Article No.	Page	Article No.	Page	Article No.	Page
JG786U	82	JK262R	46	JK487	32
JG787B	82	JK267R	46	JK488	32
JG787U	82	JK268R	46	JK489	32, 33
JG788B	82	JK286	40	JK490	44
JG788U	82	JK288	40	JK740	28 ff
JG790B	82	JK289	40	JK741	28 ff
JG790U	82	JK290	40	JK742	28 ff
JG792B	82	JK339	30 ff	JK744	28 ff
JG792U	82	JK340	30 ff	JK775A	34, 35
JG793B	82	JK341	30 ff	JK776A	34
JG793U	82	JK342	30 ff	JK777A	34
JG794B	82	JK344	30 ff	JK778A	34
JG794U	82	JK346	30 ff	JK779A	34, 35
		JK375A	36, 37	JK785	34, 35
JH217R	55	JK376A	36	JK786	34
		JK377A	36	JK787	34
JJ586B	84	JK378A	36	JK788	34
JJ586U	84	JK379A	36, 37	JK789	34, 35
JJ587B	84	JK385	36, 37		
JJ587U	84	JK386	36	JM020	43
		JK387	36	JM021	43
JK020	42	JK388	36	JM174	39
JK060R	46	JK389	36, 37	JM188	39
JK062R	46	JK395	36		
JK067R	46	JK396	36	JN021	42
JK068R	46	JK397	36	JN086	40
JK089	88	JK398	36	JN088	40
JK090	88	JK399	36, 37	JN089	40
JK091	42, 88	JK440	26 ff	JN090	40
JK092	88	JK441	26 ff	JN091	40
JK095	89	JK442	26 ff	JN092	40
JK096	86	JK443	44	JN094	40
JK098	43, 86	JK444	26 ff	JN095	40
JK100	86	JK446	26 ff	JN096	40
JK170	39	JK475A	32, 33	JN187	39
JK171	39	JK476A	32	JN188	39
JK172	39	JK477A	32	JN292	40
JK173	39	JK478A	32	JN294	40
JK174	39	JK479A	32, 33	JN295	40
JK187	39	JK485	32, 33	JN296	40
JK188	39	JK486	32	JN340	36, 37

Article No.	Page
-------------	------

JN341	36, 37
JN342	36, 37
JN344	36, 37
JN346	36, 37
JN440	32, 33
JN441	32, 33
JN442	32, 33
JN443	44
JN444	32, 33
JN445	44
JN446	32, 33
JN740	34, 35
JN741	34, 35
JN742	34, 35
JN744	34, 35

JP050	86
JP055	86
JP101	26, 27
JP102	26
JP103	26
JP104	26
JP105	26, 27
JP111	28, 29
JP112	28
JP113	28
JP114	28
JP115	28, 29
JP121	30, 31
JP122	30
JP123	30
JP124	30
JP125	30, 31
JP131	30, 31
JP132	30
JP133	30
JP134	30
JP135	30, 31

PL960R	61
PL961R	44, 61

PL963R	62
PL964R	62

AESCULAP® – a B. Braun brand

Aesculap AG | Am Aesculap-Platz | 78532 Tuttlingen | Germany
Phone +49 7461 95-0 | Fax +49 7461 95-2600 | www.aesculap.com

The main product trademark "Aesculap" and the product trademarks "PRIMELine", "instacount", "Aesculap InOS" and "STERILIT I" are registered trademarks of Aesculap AG. "Sterrads 100 S", "Sterrads 200", "Sterrads NX" and "Sterrads 100 NX" are registered trademarks of ASP. "V-Pro 1" and "V-Pro 1 Plus" are registered trademarks of Steris.

Subject to technical changes. All rights reserved. This brochure may only be used for the exclusive purpose of obtaining information about our products. Reproduction in any form partial or otherwise is not permitted.